Handleiding
Directions for use
Notice d'utilisation
Gebrauchsanweisung

Connect Carrier Cargo

CE

MD

ROAM Carrier Solutions B.V.
Buitenkamp 20
5321KB Hedel
The Netherlands
T: 0031 (0)73-5997042
E-mail: info@roam.nl
W: www.roam.nl

Inleiding/Introduction/Introduction/Einleitung

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Since 1979, Roam Carrier Solutions has been developing mobility solutions for active wheelchair users. At Roam, we offer a wide range of mobility solutions that support an active lifestyle, with durability and high quality as our guiding principles.
Every day, our team is fully committed to delivering the best mobility solutions for our customers.
#Ride that moves you!

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The management system of ROAM CARRIER SOLUTIONS BV is certified and complies with the EN ISO 13485 standard. ROAM CARRIER SOLUTIONS BV declares, as manufacturer, that this product complies with the EU Medical Device Regulation (MDR) (2017/745).

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Table of contents

Introduction	2
Table of Contents	35
Definitions	36
Delivery inspection	36-37
General	37
Applications	38
Regulations	39
Use	40-42
Attaching and detaching the Connect Carrier Cargo	42
Controls	43
Battery and charger	44-49
Display	50-53
Waste disposal	54
Reuse	55-56
Maintenance and cleaning	57
Technical specifications and characteristics	58
Adjustability	59
Possible problems	60-61
Warranty conditions	61-62
Frame sticker	62
Manufacturer contact details	63
Closing statement	63

Definitions

Always observe the definitions below!

Word	Definition
 **DANGER!**	The user is alerted to the risk of serious injury or death if the advice is not followed.
 **WARNING!**	The user is alerted to the risk of injury if the advice is not followed.
 **CAUTION!**	The user is alerted to the risk of damage to the equipment if the advice is not followed.
NOTE:	General advice or best practice
	Reference to additional documentation

Delivery check

- Connect Carrier Cargo
- Roam Adapter (if applicable)
- 2 chargers
- 2 batteries
- 1 separate rear lighting unit
- Manual
- Other components depend on the options selected.

Check whether the delivery matches the delivery note. If transport damage is found upon delivery, report this (within 24 hours) to your Roam dealer, supported with photos. Roam trusts that the Connect Carrier Cargo meets your expectations and wishes you enjoyable cycling.

NOTE!:
Before delivery of your Roam, all checkpoints have been inspected using a checklist. Certain parts must be readjusted after some time. We recommend having this carried out by your Roam dealer.

NOTE!:
Write down the details of the dealer who supplied your product here. In the event of malfunctions, contact your dealer.

Dealer name:

Dealer phone number:

Dealer email address:

Date:

General

NOTE!:
Failure to follow these instructions may result in physical injury, damage to the product, or damage to the surroundings! A notice to the user and/or patient: If this product has been involved in a serious incident, it must be reported to Roam Carrier Solutions B.V. and the competent authority of the Member State in which the user and/or patient is established.
The expected service life of the product is 7 years.
It may be possible to combine this medical device with one or more other medical devices or another product. For more information, please contact Roam Carrier Solutions B.V.

No modifications to the Connect Carrier Cargo

**WARNING!**
Roam cannot be held liable for incorrect use or incorrect information regarding the use of the Connect Carrier Cargo, and any resulting damage, unsafe situations/actions, accidents or incidents.
No modifications may be made to the Connect Carrier Cargo.
Roam is not responsible for damage to the Connect Carrier Cargo caused by modifications.
The speed of our products must not be increased or altered in any way. If the product is modified, the producer/manufacturer is in no way responsible or liable.

Applications

Indications
The Connect Carrier Cargo provides relief and compensation in cases of physical immobility, disability, or another limitation. The Connect Carrier Cargo is a Class I medical device, in accordance with Annex VIII of the MDR.

Contra-indications
The Connect Carrier Cargo must not be used without supervision in cases of:

- Perceptual disorders
- Loss of both arms
- Inability to adequately control braking and steering
- Recent surgery to the arms, shoulders and hands
- Severe fatigue or limited capacity

Additional
NOTE!:
The user must be able to operate the Connect Carrier Cargo properly.
The user of the Connect Carrier Cargo must be in good mental health and must not be physically restricted in a seated position.
The rider's hearing and eyesight must be adequate for using the Connect Carrier Cargo.
The user must read and understand the manual thoroughly.
If the user is so limited that they have vision/hearing problems and/or cannot read and understand the manual, parents/carers/supervisors are responsible for applying the manual and for the end user's use of the Connect Carrier Cargo.
The Connect Carrier Cargo is for personal use. Only the wheelchair user is permitted to use the Connect Carrier Cargo. The user is responsible and liable for correct use as stated in this manual.

It is prohibited to use the Connect Carrier Cargo under the influence of alcohol/drugs or mind-altering medication.
Read the instructions for use of the Connect Carrier Cargo. Instruct the user on the safe use of the Connect Carrier Cargo and any additionally fitted components. Inform users about specific warnings that must be read, understood and observed. When delivery/purchase of a Connect Carrier Cargo is being considered, take into account the following: body size, weight, physical and mental condition, age, living conditions and environment.

Important: Do not use your Connect Carrier Cargo until you have read and understood these instructions for use.

Regulations

NOTE!:

The user must at all times comply with the laws and regulations applicable in the country of use regarding road safety and the use of this product. The manufacturer accepts no liability for damage resulting from failure to comply with these regulations.

NOTE!:

For users under 16 years of age, it is recommended that the Connect Carrier Cargo is used under the supervision of a parent or carer.
The user may use: the footpath, the (moped) cycle path and the roadway.
It is prohibited to use motorways or expressways.
Comply with the traffic rules applicable to the road you are using. On the pavement, the maximum speed is 6 km/h.
Electronic devices such as mobile phones, navigation systems, tablets or music players must not be held while riding; Roam strongly discourages using them while riding. The same rules apply as for texting, calling or listening to music while cycling.

Use

WARNING!

1. Adjust your speed to the conditions (e.g. heavy traffic, poor road surface, obstacles and roads alongside water). Cross speed bumps and other unevenness at an appropriate speed.
2. Keep both hands on the handlebars/grips while riding and do not allow yourself to be distracted. Never put your hands into rotating parts.
3. Only shift while riding without applying force to the grips. Shift one gear at a time.
4. Be extra alert around obstacles, raised edges, railway crossings, tram tracks, zebra crossings, and when riding along busy roads or water.
5. Use lights at dusk, in the dark or in poor weather conditions. Roam recommends using lights during the day as well.
6. Roam recommends that wheelchair users wear a helmet when using the Connect Carrier Cargo.
7. Make sure you are well prepared when travelling: keep the contact details of your Roam dealer to hand (e.g. in your mobile phone).
8. On longer rides, protect yourself against cold and precipitation in winter and against excessive heat in summer.
9. Before every ride, check that the adapter and/or clamping systems under the wheelchair and/or on the Connect Carrier Cargo are fastened correctly. With detachable adapters, check the star knobs/quick releases.
10. Have maintenance, adjustments, replacement of parts and repairs carried out by your Roam dealer.
11. Take right-angled and sharp turns at an appropriate speed to prevent tipping.

Slopes and brakes

WARNING!

1. In some cases, if this manual is not followed correctly, a dangerous situation may arise for the Connect Carrier Cargo user when riding on slopes.
2. We recommend not using the Connect Carrier Cargo on slopes of 5% or more. On very long or steep slopes, the brakes may overheat and the motor may be overloaded, which can cause defects.
3. The point above regarding slopes also applies to descents.
4. Counterweights (optional) can be fitted, generating more wheel load and therefore more grip on the front wheel. This advantage of increased wheel load is beneficial for uphill slopes, not for descents (speed may increase).
5. Be aware that speed increases when descending a slope.
6. Please note: On roads with oil, mud, sand or leaves on the surface, the braking distance may be longer. This also applies to uneven road surfaces (potholes/holes) and winter conditions.

Corners

DANGER!

Take right-angled and sharp turns at an appropriate speed, depending on the conditions (road surface and weather) and your physical capabilities, to prevent tipping. Make sure the grips are in the upper position before taking the turn so you can steer easily. This prevents your thighs or upper body from being hit by the grips, which could restrict the steering angle.

Maximum user weight

DANGER!

Never exceed the maximum user weight of 170 kg. This is separate from the maximum user weight of the wheelchair; please consult the relevant wheelchair supplier/manufacturer for that information. Exceeding the maximum load may result in damage to the Connect Carrier Cargo, or it may cause tipping or overturning and loss of control, which could lead to serious injury to the user and other persons.

Choking hazard

WARNING!

This mobility aid contains small parts which, under certain circumstances, may pose a choking hazard for small children. The product shown and described in this manual may not be exactly the same in every detail as your model. However, all instructions are important regardless of differences in detail. The manufacturer reserves the right to make changes without prior notice to information regarding weight, dimensions or other technical data as stated in this manual.

Wheelchair

WARNING!

If the wheelchair's balance point is set very sharply/tip-sensitive, the Connect Carrier Cargo will have reduced wheel load and the braking distance will be much longer. Reduced wheel load can cause the risk of slipping/wheel spin, including uphill.

NOTE!:

The Roam adapter causes the wheelchair to become 1 to 2 kg heavier. This depends on the type of wheelchair and the associated mounting components of the Roam adapter. This has little to no effect on the driving characteristics of the wheelchair itself when the Connect Carrier Cargo is not in use.
If the user wishes to remove the adapter from the wheelchair, this can be done by loosening the M6 star knobs (counterclockwise) and removing the adapter from the wheelchair with a forward motion. This applies to a detachable adapter.

The wheelchair can be fully restored to its original condition after removal of the Roam Connect Carrier Cargo and the associated adapter assembly. All fasteners and mounting components must be completely removed from the wheelchair frame.

Attaching and detaching the Connect Carrier Cargo

Attaching using the adapter system
Connecting the Connect Carrier Cargo to your wheelchair is very simple. First, apply the wheelchair brake; then close the connector lever on the Connect Carrier Cargo (downward) and slide the Connect Carrier Cargo with the coupling frame as far as possible into the adapter on your wheelchair. Now open the connector lever on the Connect Carrier Cargo (upward), so that the bracket falls over the hook section of the adapter, then close the lever again. The Connect Carrier Cargo is now connected. Release the wheelchair brake and you can start driving.

NOTE!:
The actions for connecting and disconnecting may vary slightly depending on the options selected. This manual is based on the most basic version. Please refer to the instruction videos on the website www.roam.nl

Detaching using the adapter system
Disconnecting the Connect Carrier Cargo is also very simple. First, apply the wheelchair brake, then hold the Connect Carrier Cargo, preferably by gripping the handlebar. Next, use your other hand to open the connector and hold the clamp bracket upward. You can now close the connector lever again so that the bracket rests on the hook section. The Connect Carrier Cargo can then be slowly slid out of the adapter in order to set it down.

NOTE!:
The actions for connecting and disconnecting may vary slightly depending on any selected options. This manual is based on the standard configuration in the most basic version. Please refer to the instruction videos on the website www.roam.nl

Controls

Brakes

DANGER!

It is important that you fully understand how the braking system of your Connect Carrier Cargo works. If the braking system is not used correctly, you may lose control or it may lead to serious consequences (accident, serious injury).
The Connect Carrier Cargo is equipped with a disc brake. If, in the user's opinion, the brake does not function properly, contact the dealer immediately. The Roam dealer must resolve any defects. It is prohibited to use the Connect Carrier Cargo in any way if the braking system is not functioning.

WARNING!

1. The disc brake can become hot; contact the Roam dealer immediately if there is doubt or improper brake function.
2. Do not come into contact with the disc brake, frame parts or housing.
3. Always check your brake before riding.
4. It is prohibited to use oil or grease near the disc brake system.

Parking brake function

When attaching and detaching, the Connect Carrier Cargo must not be on the handbrake. After detaching, the Connect Carrier Cargo must always be put back on the handbrake. Always use the parking brake when loading and unloading.
Pull the large lever towards you and then squeeze the brake lever; you will hear a "click" sound. Do this until you can no longer squeeze the brake lever any further. The Connect Carrier Cargo is now on the parking brake.
To release the parking brake, swing the lever on the brake lever away from you; the parking brake is then released.

WARNING!

If the slope gradient is more than 5%, the brake can become hot and may no longer brake effectively! In this situation you must keep your speed limited (no more than 20 km per hour).

Speed
NOTE!:
The Connect Carrier is supplied as standard with 5 speed settings, which can be selected via the display. The maximum speed of the Connect Carrier Cargo is 18 km/h.

CAUTION!

- The propulsion speed depends on the user, the surface, and the operating conditions.
- Use the aid only at a safe and controlled speed.
- Avoid high speed on slopes and on wet or uneven surfaces.
- The maximum safe speed is determined by the user's control and the environment.
- Sudden changes in direction at higher speed may lead to loss of control.

Lighting

CAUTION!
The Connect Carrier Cargo is equipped as standard with front lighting, rear lighting, and a separate rear lighting unit.
The rear light must be mounted on the back of the relevant wheelchair frame using the supplied plastic bracket with elastic strap mounting.
Switch on the rear light by pressing the button on the rear lighting unit.

Battery and charger

WARNING!
Read and follow the safety instructions before using the battery and before starting the charging process.

DANGER!
The lithium-ion battery contains chemicals that may cause hazardous reactions if the safety instructions specified in this manual are ignored.

NOTE:
Roam Carrier Solutions accepts no liability for damage resulting from failure to comply with these instructions.

Safety instructions and precautions for battery use

1. Before using the Connect Carrier Cargo for the first time, the batteries must be fully charged.
2. The batteries may only be used between -15°C and 45°C.

3. The batteries must not be exposed to heat or open flame. External heat may cause the batteries to explode.
4. In the (unlikely) event that the battery overheats or catches fire, the battery must absolutely NOT come into contact with water or other liquids.
5. The batteries must not be opened or dismantled. Unauthorized opening or deliberate damage to the batteries carries a risk of serious injury. In addition, any warranty claim becomes void if the battery is opened.
6. Never connect the battery terminals in the charging socket to metal objects, and ensure that the terminals never come into contact with metal objects.
7. The service life of the battery depends in part on where it is stored. Never leave the battery in a hot location for an extended period. In particular, the luggage compartment of a car parked in the sun can be hazardous. Use the luggage compartment only for transport and not as a storage location for the battery.
8. The battery must not be subjected to mechanical impact. If this does occur, the battery must be inspected by the manufacturer. Please contact your dealer about this. You must not continue driving with a damaged battery.
9. Always ensure that the batteries are stored dry and clean.
10. Place the batteries on a fire-resistant and heat-resistant surface during charging. Do not place easily flammable or combustible objects close to the battery.
11. During charging, do not place the charger on a wet surface or outside in the rain.
12. Do not use the batteries or chargers if plugs or connectors are wet.
13. Improper use may lead to overheating, fire, or explosion.
14. The battery must be charged to at least 20% before use.
15. Ensure that the battery is not discharged below 20% during use.
16. Store the battery fully charged; if it is not used for an extended period, repeat this every 3 months.
17. The ideal storage temperature for the battery is between +10°C and +25°C, with a humidity level between 20 and 75 percent. This ensures the optimum service life of the battery.
18. Use only original Roam chargers to charge the battery.
19. Charge the batteries in a safe place under supervision; never charge them at night due to the fire hazard in the event of overheating.
20. Replace damaged chargers through your Roam dealer.
21. Do not store batteries in direct sunlight or behind a window in the sun.
22. When the charging cycle is complete, first unplug the charger from the power outlet before removing the battery from the charger.
23. Ensure adequate ventilation when the batteries are being charged.

NOTE:*We recommend charging the batteries in the fire-resistant case, which is available as an optional accessory. In addition, we recommend storing the batteries in the fire-resistant case.

Transport
If your battery is defective, take it personally to your specialist dealer; sending lithium-ion batteries by post or other means of transport is subject to strict laws and regulations. Once again, we advise you to contact your dealer.

As transport regulations change every year, we strongly recommend consulting your travel agency, airline or shipping company before travelling, so that you know the current applicable regulations. A defective battery must not be taken on an aircraft or placed in your luggage.

NOTE!: Discuss transport with your dealer prior to shipment.

Using the Roam Smart charger

The MDA BC257360020 is a Li-ion battery charger with:

- Input: 100–240V AC
- Output: approx. 42V DC, 2A
- Application: charging 36V lithium-ion batteries

Use (step-by-step)

1. Check the battery
 - Only charge if it is compatible with the charger (Li-ion).
 - Do not use a damaged battery.
2. Connect the battery first
 - Connect the charger plug to the battery.
3. Then plug the charger into the mains socket
4. Indicator lights
 - Red = charging in progress
 - Green = fully charged (may vary slightly depending on the version)
5. After charging
 - First unplug from the mains socket
 - Then disconnect from the battery

Important safety instructions

- For indoor use only (no rain or moisture)
- Do not leave the charger unattended during use
- Do not open or attempt to repair yourself
- Use only with compatible batteries
- Keep out of reach of children
- Do not use in case of:
 - smoke
 - unusual smell
 - damaged cable
- Unplug from the mains socket before cleaning

Protective functions (automatic)

This charger typically has built-in protections such as:

- Overcharge protection
- Short-circuit protection
- Overheating protection
- Overvoltage protection

Tips for optimal use

- Fully charge the battery before first use
- Use the battery regularly (prevents loss of capacity)
- Do not store batteries empty
- Preferably charge at room temperature

Common problems

The charger does not work

- Check the mains socket
- Check the cables

The battery does not charge

- Possibly incorrect battery
- Defective or deeply discharged battery

The charger becomes warm

- Slightly warm = normal
- Very hot = stop immediately

Battery 36V/14.4Ah/518.4Wh

This manual describes:

- Product features and performance
- Correct use of the battery
- Safety instructions for lithium-ion battery packs

Product overview

Main parts

The battery consists of the following components:

- Lock incl. 3 keys
- LED indicator
- On/off button
- Discharge connector (to the bike)
- Locking lever (drop protection)
- Charging connector
- Warranty seal
- Product label

Technical specifications

- Battery type: Lithium-ion
- Nominal voltage: 36V
- Capacity: 14.4Ah

Use

Switching the battery on and off

- Press the on/off button to activate the battery
- The LED indicators show the charge level

Installing the battery

- Place the battery in the bike's holder
- Ensure it clicks firmly into place

Charging

First connect the charger to the battery.

1. Then plug the charger into the wall socket
2. Charging starts automatically
3. After fully charging:
 - First unplug the charger from the wall socket
 - Then disconnect the battery

LED indicator

- Multiple LEDs lit = charge level
- Flashing LEDs = charging in progress
- All LEDs lit = fully charged
 (The exact pattern may vary slightly depending on the version.)

Safety instructions

General safety

- Use only the specified charger
- Do not disassemble or modify the battery
- Protect against water and fire
- Avoid impacts, drops and damage

Safe charging

- Charge indoors only, in a dry place
- Do not cover the battery while charging
- Stop using immediately if you notice:
 - unusual smell
 - overheating
 - deformation

Safe use

- Prevent short circuits
- Do not expose to temperatures above 60°C
- Keep out of reach of children

Storage

- Store in a cool, dry place
- Ideal state of charge: 40–60%
- Recharge at least every 3 months
- Do not store empty

Maintenance

- Keep connectors clean and dry
- Do not use aggressive cleaning agents
- Check regularly for damage

Common problems

Battery does not charge

- Check the charger and connection
- Check the wall socket
- The battery may be deeply discharged

Battery does not work on the bike

- Check that it is switched on
- Check that it is installed correctly

Battery becomes too hot

- Stop using it immediately
- Allow the battery to cool down

Disposal

- Do not dispose of with household waste
- Follow local regulations for lithium batteries
- Take it to an authorised collection point

Display

CAUTION!

Do not connect or disconnect the display while it is powered. This may damage electronic components.

Due to product updates, the display may differ slightly from this manual, without affecting normal use.

Technical specifications

- Power supply: DC 24V / 36V / 48V
- Rated current: 23 mA
- Leakage current (off): < 1 μA
- Display: 2.0" IPS LCD, resolution 320 × 480
- Communication: UART / CAN bus
- Operating temperature: -20 °C to 60 °C
- Storage temperature: -30 °C to 80 °C
- Waterproof rating: IP65

Key functions

- Ergonomic buttons on the left side
- Customizable start-up and user interface
- Units: km/mile, languages: English/German
- Display of ride data (speed, distance, battery, etc.)
- Assist levels (power assist) with statistics
- Walk assist function
- Automatic lighting (on/off)
- Error messages with codes
- Built-in clock (RTC)
- Range and battery indication (depending on BMS)
- Power distribution between rider and motor (if supported)
- Maintenance reminder
- Advanced settings (for professionals)

Installation (brief)

- Secure with M3 screw (max. 0.8 Nm)
- Connect via 5-pin connector to the controller
- Note: overtightening may cause damage

Use & operation

Start-up

- Display shows the start-up screen for 3 seconds
- Then the main interface appears with real-time data

Buttons

- On/off button: switch the display on/off
- M button: switch between screens / settings
- + / – buttons: adjust assist level

Basic operation

- Short press = activate function
- Long press (±2 sec) = additional function

Main screens

- Basic interface: speed, battery, assist level
- Screen 1: average/max speed and trip distance
- Screen 2: battery data and range
- Screen 3: motor power and power distribution
- Screen 4: assist level usage statistics

Use the M button to switch between screens

After 5 sec without action, the screen returns to the basic interface

Important functions

Assist level

- Adjustable from 0 to 5
- Not cyclic (no automatic reset to 0)

Walk assist

- Press and hold when stationary
- Motor provides support for slow walking

Lighting

- Automatic via light sensor
- Manual on/off via button
- Screen dims when the light is switched on

Maintenance reminder

- Standard after 5000 km
- Can be turned on/off in settings

Battery display

- 5 levels + percentage (if available via BMS)
- Below <5%: warning + no motor assistance
- In case of communication problems:
 - Estimate only based on voltage
 - No percentage or range displayed

Settings

- Within 10 sec after start-up: press and hold M
- Navigation:

Short M = select
/ – = adjust
Long M = confirm / back

Advanced settings

- For manufacturers/dealers only
- Password required
- Incorrect use may cause system malfunctions

Reset data

- Reset possible for:
 - Trip distance
 - Average and maximum speed


Odometer (total km) cannot be reset manually

Error codes

Error code	Error description	Suggested action
"04" displayed on speed	Throttle does not return to zero position (remains in the high position)	Check whether the throttle has returned
"05" displayed on speed	Throttle faulty	Check throttle
"07" displayed on speed	Overvoltage protection	Check the battery voltage
"08" displayed on speed	Motor Hall signal wire fault	Check motor
"09" displayed on speed	Motor phase wire fault	Check motor
"11" displayed on speed	Motor temperature sensor fault	Check controller
"12" displayed on speed	Current sensor fault	Check controller
"13" displayed on speed	Battery temperature failure	Check battery
"14" weergegeven op snelheid	The controller temperature is too high and reaches the protection point	Check motor
"21" displayed on speed	Speed sensor failure	Check the installation position of the speed sensor
"22" displayed on speed	BMS communication fault	Replace battery
"30" displayed on speed	Communication error	Check connector to controller

Waste disposal

Below are the materials as used in the Connect Carrier Cargo, with regard to waste disposal and recycling of the Connect Carrier Cargo and the packaging.
Specific applicable legislation regarding waste disposal and recycling must be observed when the Connect Carrier Cargo has reached the end of its service life.
Aluminium: steering tube, rim
Steel: front fork, stem, motor
Stainless steel: stand, coupling frame
Plastic: grips, chain guard, (inner) tube
Packaging: cardboard, plastic

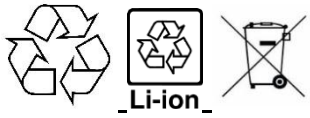Reuse

The Connect Carrier Cargo is suitable for multiple use, provided it is cleaned and complies with the inspection procedures.
Check the following components for wear, damage, and functionality before reuse:

- Batteries: test these. If the measured battery capacity is above 70%, the battery is approved; below 70%, it is rejected. Roam can perform battery testing.
- Inspect the frame for visible cracks and paint damage.
- Inspect all brake cables for damage and replace them if necessary.
- Inspect the wheel for loose spokes and play in the wheel bearings.
- Inspect the tyre for wear and replace it if necessary.
- Inspect the electrical cables and plug connections for cable breakage, visible damage, and corrosion, and replace them if necessary.
- Inspect all 4 brake calipers and the 2 brake discs for proper function, including the corresponding brake pads. Replace the entire brake caliper, brake discs, or brake pads in the event of poor performance or visible wear.
- Inspect the gears for proper function.
- Inspect which options have been added to the Connect Carrier Cargo and adjust or replace them if necessary.
- Adjust the Roam adapter where necessary so that it fits the other wheelchair. NOTE: This must at all times be carried out by your mobility aid dealer.
- Inspect the connector for preload and adjust it if necessary; see the instruction sheet below for the connector and/or connector lever assembly instructions. If components show visible wear, always replace the complete connector lever assembly (complete), item 100318, or connector (complete), item no. 9702022.


Instruction sheet


1. Loosen the M6 bolts and slide over the coupling tube; if necessary, the connector lever can be removed. Slide to the correct position and mount straight relative to the wheelchair. Tightening torque for M6: 9-11 Nm

2. Adjust the connector lever by setting the M6 nuts.

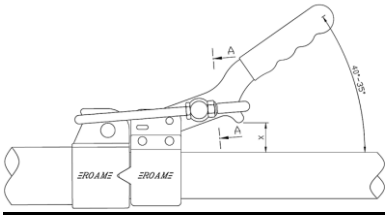3. The connector is correctly adjusted if:
 - the lever starts clamping from 40 - 35º relative to the coupling tube and/or
 - the distance X = is between 20 and 15 mm

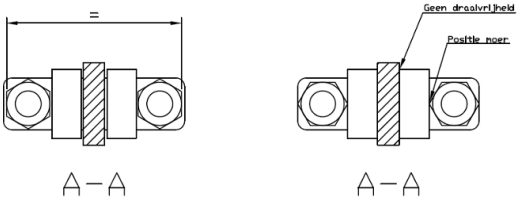[section A-A] Ensure that the M6 nuts are flush on the flat side. This ensures rotational clearance at the plastic bush, allowing the lever to move smoothly.

Maintenance and cleaning

CAUTION!

Vital parts such as cables, spokes and brakes may settle during the initial period. To keep the Connect Carrier Cargo in good condition, an annual service is recommended (necessary). Your Roam dealer has the correct tools and experience to adjust and lubricate hidden rotating parts such as the headset and hubs. Good maintenance is also extremely important with regard to the warranty. Make an appointment with your Roam dealer for servicing.

Cleaning

1. Clean the Connect Carrier Cargo with a dry cloth.
2. If very dirty, clean the Connect Carrier Cargo with a soft, damp sponge and dry afterwards.
3. Do not spray the Runner Connect Carrier Cargo with a high-pressure washer.
4. Do not use chemical cleaning agents (chlorine/bleach/soda) when cleaning the Connect Carrier Cargo.
5. Preferably store the Connect Carrier Cargo dry or at least under cover to prevent oxidation.

Tools intended for Roam dealer

The mechanic at your Roam dealer has the correct tools to maintain/repair the Connect Carrier Cargo. Original Roam parts must be used for replacement. Your Roam dealer can order these from Roam in Hedel. Ask your Roam dealer for advice.

NOTE!:

Check regularly before using the Connect Carrier Cargo:

1. Check brake function.
2. Check wear of the brake pads.
3. Check the tightening torque of all bolts and nuts.
4. Check the condition of the tyre for bumps, cracks and tread depth.
5. Check tyre pressure (recommended 2.5–3.5 bar).
6. Check the function of the front and rear lights.
7. Check the sidewall reflectors of the outer tyre.
8. Check that the grips are mounted correctly and still provide the correct grip.
9. Check the operation of the steering corrector (if fitted).
10. Check the clamping force of the connector.
11. Check the mounting of the adapter under the wheelchair frame.

WARNING!

If you have any doubts about the operation of the Connect Carrier Cargo, contact your Roam dealer immediately.

Technical data and characteristics

Adjustable coupling tube:	stainless steel
Connect Carrier Cargo frame:	steel
Stem:	Roam adjustable aluminium stem
Rim:	Ryde Andra 29 - 20 inch
Spokes:	Saphim
Motor:	Roam
Band:	Schwalbe Pick-Up 55-406 20'' x 2.15
Brakes:	Tektro Auriga Twin HD 545
Front light:	Spanninga Pimento led (StVZO approved)
Rear light:	Spanninga Pimento Brake led (StVZO approved)
Motor:	2x 36V - 250W
Battery:	2x 14,4Ah. / 2x 520Wh. / 36V Lithium Ion
Overall width:	83 cm
Overall height:	86 cm
Overall length:	185 cm incl. coupling frame
Weight in kg:	52 kg excl. battery
Maximum load:	170 kg.
Maximum payload:	50 kg.
Cargo space volume:	195 litres
Wheel size:	20 inch
Maximum speed:	18 km/h

Tightening torques
Below you will find an overview of the recommended tightening torques. We recommend tightening the listed components with a torque wrench according to the specified values. Overtightening the components may cause damage to your Connect Carrier Cargo. For this procedure, it is best to ask your Roam dealer for help and advice.

- Expander bolt with tapered wedge 12–15 Nm
- Front wheel hub axle nuts 20–27 Nm
- Connector bolts (connector at wheelchair adapter) 4–5 Nm
- Connector bolts (connector at Connect Carrier Cargo coupling tube) 9–11 Nm
- Connector locking bolt 9–11 Nm (apply thread-locking compound)
- Coupling frame/head tube aluminium clamp 38 mm 9–11 Nm
- Plastic adapter clamp block bolts 9–11 Nm

NOTE!:
When loosening the M6 x 20 mm hex socket bolts of the Roam connector on the coupling frame after making an adjustment, it is recommended to use new M6 x 20 mm bolts and retighten them with assembly grease using a torque wrench to 9–11 Nm, and tighten the connector clamped at the wheelchair adapter to 4–5 Nm.

Adjustability


The following components can be adjusted to change the user's seating position or to adjust the settings of the Connect Carrier Cargo in combination with the wheelchair.

Component:	Adjustment direction:	Effect:
Roam adjustable stem	Height and length	Handlebar position
Connector coupling frame	Lengthwise	Carrier position and front wheel lift
Connector preload	Connector clamping tension	Can be adjusted tighter or looser

Possible problems

- Connect Carrier Cargo is misaligned in front of the wheelchair:
 Loosen the 2x M6 bolts of the connector (cast part) on the coupling frame and connect the Connect Carrier Cargo to the wheelchair adapter. Align the Connect Carrier Cargo straight and retighten these bolts alternately.
- Check whether the adapter is mounted in the centre of the wheelchair and whether the bolts are tightened to the correct torque.
- Squeaking or scraping noise at the wheel:
 A dragging brake caliper may be the cause; readjust the brake caliper and brake cable.
 The brake pads are worn; in that case, replace them immediately.
- Play in the connector:
 Readjust the connector lever. Ask the relevant dealer to have this carried out.
- Connect Carrier Cargo does not connect properly in combination with the adapter:
 Check whether the adapter is mounted correctly under the wheelchair. It may have shifted, causing the Connect Carrier Cargo not to couple and uncouple properly.
 The wheelchair settings have been changed, causing the Connect Carrier Cargo not to couple and uncouple properly.

NOTE!:
The connector (see image) must be able to move freely widthwise over the coupling tube.
The connector is secured on the underside with 1 pc. M6x20 mm (see image) to prevent it from moving too far. This M6 bolt fits into a hole in the stainless steel tube; this hole has a diameter of approximately 8 mm.
Tighten the bolts according to the tightening torques stated in the drawing below.

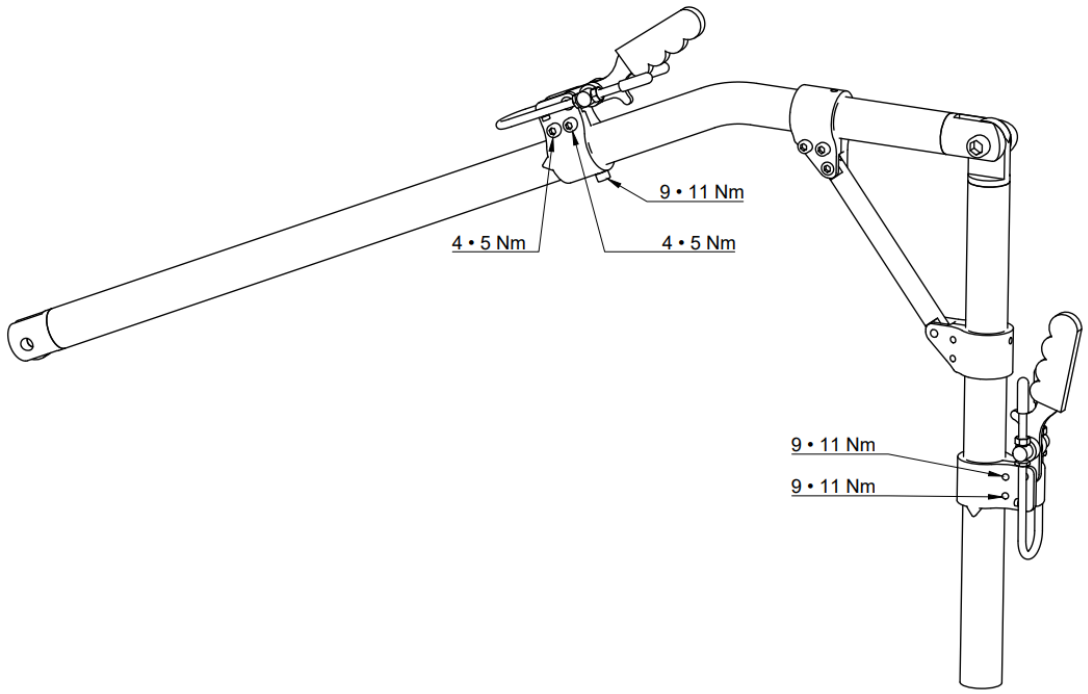Warranty Terms and Conditions

Roam Special Cycles BV provides a lifetime warranty against manufacturing, material and construction defects for the frame and front fork, calculated from the date of purchase.

All parts and components manufactured and marketed by Roam Special Cycles BV are covered by a 2-year warranty against material and manufacturing defects, calculated from the date of purchase by the end user, unless stated otherwise.

The warranty excludes:
• normal wear and tear;
• tyres;
• wear parts such as brake pads, brake discs, chains, sprockets, bearings and other parts subject to normal use;
• damage resulting from improper use, insufficient maintenance, accidents, modifications to the product or use contrary to the user manual.

If a defect covered by the warranty occurs within the applicable warranty period, Roam Special Cycles BV will, at its own discretion, repair or replace the relevant part. Any shipping or transport

costs for warranty parts will not be charged, unless parts unrelated to the warranty claim are shipped in the same consignment.


Warranty claims from end users must be submitted through an official Roam dealer. The relevant Roam dealer registers the product with Roam Special Cycles BV by email, stating a clear description of the complaint, the packing slip or invoice number and the relevant part. After registration, the Roam dealer sends the relevant part, together with a fully completed RMA form, to Roam Special Cycles BV.

After receipt by Roam Special Cycles BV, the warranty claim will be assessed and processed within 3 to 5 working days, unless communicated otherwise.

Please note: a copy of the original proof of purchase or packing slip and a fully completed RMA form must always be submitted. This form can be downloaded from www.roam.nl.

The warranty applies exclusively to the first owner of the product and is not transferable, unless agreed otherwise in writing.

Frame sticker

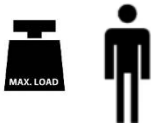TYPE:	Product name/SKU number.
	Maximum user weight.
CE	CE conformity marking.
MD	Medical device.
UDI	Unique device identification.
SN	Serial number.
	Name and address of the manufacturer.
 YYYY-MM-DD	Production date.
	Consult the instructions for use.

Manufacturer contact details

Roam Carrier Solutions B.V.
Buitenkamp 20
5321 KB Hedel
T: 0031 (0)73-5997042
E-mail: info@roam.nl
W: www.roam.nl

Closing statement

Thank you for choosing a product from Roam Carrier Solutions BV. This product has been developed in accordance with high quality and safety requirements. With correct use and maintenance, it will contribute to reliable and comfortable mobility.
For additional information or support, please first contact your supplier/dealer or Roam.

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