Handleiding
Directions for use
Notice d'utilisation
Gebrauchsanweisung

Runner Handbike Kids Elektro 16''

CE

MD

ROAM Special Cycles B.V.
Buitenkamp 20
5321KB Hedel
The Netherlands
T: 0031 (0)73-5997042
E-mail: info@roam.nl
W: www.roam.nl

Inleiding/Introduction/Introduction/Einleitung

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Table of Contents

Introduction	2
Table of Contents	38
Definitions	39
Delivery inspection	39-40
General	40
Applications	41
Regulations	42
Use	43-46
Coupling and uncoupling the handbike	46-47
Controls	47-48
Battery and charger	49-54
Display	54-57
Accessories	57
Waste disposal	57
Reuse	57-59
Maintenance and cleaning	60
Technical specifications and features	61-62
Possible problems	62
Warranty terms and conditions	63
Frame sticker	64
Manufacturer contact details	65
Closing remarks	65
Appendices	126-134

Definitions

Always observe the definitions below!

Word	Definition
 **DANGER!**	The user is alerted to the risk of serious injury or death if the advice is not followed.
 **WARNING!**	The user is alerted to the risk of injury if the advice is not followed.
 **CAUTION!**	The user is alerted to the risk of damage to the equipment if the advice is not followed.
NOTE:	General advice or best practice
	Reference to additional documentation

Delivery check

Check whether the delivery matches the delivery note. If transport damage is found upon delivery, report this (within 24 hours) to your Roam dealer, supported with photos. Roam trusts that the Runner Handbike Elektro 16'' (hereinafter referred to in this manual as: handbike) meets your expectations and wishes you enjoyable cycling.

Scope of Delivery

- Runner Elektro Handbike 16''
- Roam Adapter (if applicable)
- Charger
- Battery (check the battery type on the packing slip and on the label of the battery housing)
- Separate rear light
- Manual
- Other components depend on the options selected.

NOTE!:
Before delivery of your Roam, all checkpoints have been inspected using a checklist. Certain parts must be readjusted after some time. We recommend having this carried out by your Roam dealer.

NOTE!:
Write down the details of the dealer who supplied your product here. In the event of malfunctions, contact your dealer.

Dealer name:

Dealer phone number:

Dealer email address:

Date:

General

NOTE!:
Failure to follow these instructions may result in physical injury, damage to the product, or damage to the surroundings! A notice to the user and/or patient: If this product has been involved in a serious incident, it must be reported to Roam Special Cycles B.V. and the competent authority of the Member State in which the user and/or patient is established.
The expected service life of the product is 7 years.
It may be possible to combine this medical device with one or more other medical devices or another product. For more information, please contact **Roam Special Cycles B.V.**

No modifications to the handbike

WARNING!
Roam cannot be held liable for incorrect use or incorrect information regarding the use of the handbike, and any resulting damage, unsafe situations/actions, accidents or incidents.
No modifications may be made to the handbike.
Roam is not responsible for damage to the handbike caused by modifications.
The speed of our products must not be increased or altered in any way. If the product is modified, the producer/manufacturer is in no way responsible or liable.

Applications

Indications
The handbike provides relief and compensation in cases of physical immobility, disability, or another limitation. The Runner Handbike Elektro 16'' is a Class I medical device, in accordance with Annex VIII of the MDR. Suitable for both adults and children with a body height from 1.05 metres to 1.85 metres and an age of 5 years and older.

Contra-indications
The handbike must not be used without supervision in cases of:

- Perceptual disorders
- Loss of both arms
- Inability to adequately control braking and steering
- Recent surgery to the arms, shoulders and hands
- Severe fatigue or limited capacity

Additional
NOTE!:
The user must be able to operate the handbike properly.
The user of the handbike must be in good mental health and must not be physically restricted in a seated position.
The rider's hearing and eyesight must be adequate for using the handbike.
The user must read and understand the manual thoroughly.
If the user is so limited that they have vision/hearing problems and/or cannot read and understand the manual, parents/carers/supervisors are responsible for applying the manual and for the end user's use of the handbike.
The handbike is for personal use. Only the wheelchair user is permitted to use the handbike. The user is responsible and liable for correct use as stated in this manual.

Additional information for the dealer
NOTE!:
It is prohibited to use the handbike under the influence of alcohol/drugs or mind-altering medication. Read the instructions for use of the handbike. Instruct the user on the safe use of the handbike and any additionally fitted components. Inform users about specific warnings that must be read, understood and observed. When delivery/purchase of a handbike is being considered, take into account the following: body size, weight, physical and mental condition, age, living conditions and environment.

Important: Do not use your handbike until you have read and understood these instructions for use.

Regulations

NOTE!:

The user must at all times comply with the laws and regulations applicable in the country of use regarding road safety and the use of this product. The manufacturer accepts no liability for damage resulting from failure to comply with these regulations.

NOTE!:

For users under 16 years of age, it is recommended that the handbike is used under the supervision of a parent or carer.

The user may use: the footpath, the (moped) cycle path and the roadway.

It is prohibited to use motorways or expressways.

Comply with the traffic rules applicable to the road you are using. On the pavement, the maximum speed is 6 km/h.

Electronic devices such as mobile phones, navigation systems, tablets or music players must not be held while riding; Roam strongly discourages using them while riding. The same rules apply as for texting, calling or listening to music while cycling.

Use

WARNING!

1. Adjust your speed to the conditions (e.g. heavy traffic, poor road surface, obstacles and roads alongside water). Cross speed bumps and other unevenness at an appropriate speed.
2. Keep both hands on the handlebars/grips while riding and do not allow yourself to be distracted. Never put your hands into rotating parts.
3. Only shift while riding without applying force to the grips. Shift one gear at a time.
4. Be extra alert around obstacles, raised edges, railway crossings, tram tracks, zebra crossings, and when riding along busy roads or water.
5. Use lights at dusk, in the dark or in poor weather conditions. Roam recommends using lights during the day as well.
6. Roam recommends that wheelchair users wear a helmet when using the handbike.
7. Make sure you are well prepared when travelling: keep the contact details of your Roam dealer to hand (e.g. in your mobile phone).
8. On longer rides, protect yourself against cold and precipitation in winter and against excessive heat in summer.
9. Before every ride, check that the adapter and/or clamping systems under the wheelchair and/or on the handbike are fastened correctly. With detachable adapters, check the star knobs/quick releases.
10. Have maintenance, adjustments, replacement of parts and repairs carried out by your Roam dealer.
11. Take right-angled and sharp turns at an appropriate speed to prevent tipping.

Slopes and brakes

WARNING!


1. In some cases, if this manual is not followed correctly, a dangerous situation may arise for the handbike user when riding on slopes.
2. We recommend not using the handbike on slopes of 5% or more. On very long or steep slopes, the brakes may overheat and the motor may be overloaded, which can cause defects.
3. The point above regarding slopes also applies to descents.
4. Counterweights (optional) can be fitted, generating more wheel load and therefore more grip on the front wheel. This advantage of increased wheel load is beneficial for uphill slopes, not for descents (speed may increase).
5. Be aware that speed increases when descending a slope.
6. Please note: On roads with oil, mud, sand or leaves on the surface, the braking distance may be longer. This also applies to uneven road surfaces (potholes/holes) and winter conditions.

Front wheel traction loss

DANGER!

When using the handbike, the front wheel may slip under certain conditions, particularly when riding on slopes or on slippery surfaces. This phenomenon can affect riding performance and steering control. Although the risk of slipping can be reduced through correct use and adjustment, it cannot be completely eliminated.

- Ensure that the tyre is not overinflated (2 to 2.5 bar)
- Ride in the highest assistance mode (mode 5)
- Maintain momentum. Build up speed in advance, not only once you are already on the slope
- Where possible, shift your body slightly forward
- Install counterweights; these reduce slipping
- If necessary, fit an additional set of counterweight plates to the counterweights
- If necessary, install a cargo carrier so that carried items are mounted on the bike and not on the back of the wheelchair
- Ensure that the lift of the wheelchair's front wheels is no more than 5 cm
- Use a more active handbike adjustment, with the unit positioned further forward
- Install a second rear axle on the wheelchair (not possible on every wheelchair; consult your dealer/wheelchair manufacturer)
- Keep your speed up until you reach the top, and only then slow down.
- Take extra care on wet, loose or slippery surfaces

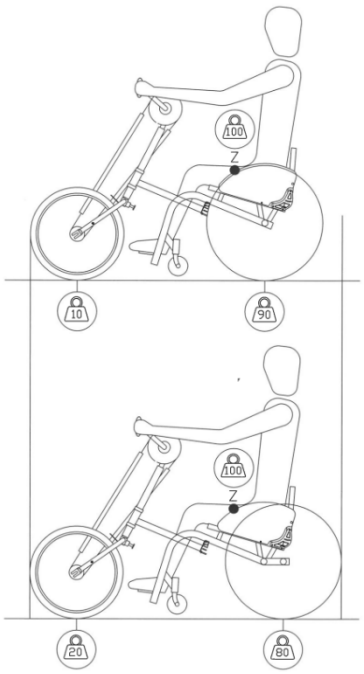*The image shows an example of a second rear axle; the pressure on the front wheel increases.

Corners

DANGER!

Take right-angled and sharp turns at an appropriate speed, depending on the conditions (road surface and weather) and your physical capabilities, to prevent tipping. Make sure the grips are in the upper position before taking the turn so you can steer easily. This prevents your thighs or upper body from being hit by the grips, which could restrict the steering angle.

Maximum user weight

DANGER!

Never exceed the maximum user weight of 170 kg. This is separate from the maximum user weight of the wheelchair; please consult the relevant wheelchair supplier/manufacturer. Exceeding the maximum load may cause damage to the handbike, or you may fall over or tip, lose control, and cause serious injury to the user and other people.

Choking hazard

WARNING!

This mobility aid contains small parts which, under certain circumstances, may pose a choking hazard for small children. The product shown and described in this manual may not be exactly the same in every detail as your model. However, all instructions are important regardless of differences in detail. The manufacturer reserves the right to make changes without prior notice to information regarding weight, dimensions or other technical data as stated in this manual.

Wheelchair

WARNING!

If the wheelchair balance point is adjusted to be sharp or tip-sensitive, the handbike will have reduced wheel load and the braking distance will be much longer. Reduced wheel load can cause a risk of slipping or wheel spin, including when riding uphill.

NOTE:

The Roam adapter increases the weight of the wheelchair by 1 to 2 kg. This depends on the type of wheelchair and the associated fastening components of the Roam adapter. This has little to no effect on the driving characteristics of the wheelchair itself when the handbike is not in use.

If the user wishes to remove the adapter from the wheelchair, this can be done by loosening the M6 star knobs (counterclockwise) and removing the adapter from the wheelchair with a forward movement. This applies to a removable adapter.
In the Klick-Fix (front frame mounting) version, nothing additional is attached to the wheelchair for coupling.

The wheelchair can be fully restored to its original condition after removal of the Roam Handbike and the associated adapter construction. In doing so, all fastening materials and mounting components must be completely removed from the wheelchair frame.
In a version with the Klick-Fix system, no permanent components or fasteners are mounted on the original wheelchair frame. As a result, the original configuration of the wheelchair remains unchanged.

Coupling and uncoupling the handbike

Coupling using the adapter system
Coupling the handbike to your wheelchair is very easy. First apply the wheelchair brake; then close the connector lever on the handbike (downwards) and slide the handbike with the coupling frame as far as possible into the adapter on your wheelchair.
Now open the connector lever on the Runner Handbike (upwards) so that the bracket falls over the tab of the adapter; then close the lever again. The handbike is now coupled. Release the wheelchair brake and you are ready to ride.

NOTE!:
The coupling/uncoupling steps may differ slightly, depending on the type of stand you have chosen. This manual assumes the bicycle stand in the most basic version. See the instruction videos on the website www.roam.nl

Uncoupling using the adapter system
Uncoupling the Runner Handbike is also very easy. First apply the wheelchair brake. Then hold the Runner Handbike (preferably by holding the handlebars); with your other hand open the connector and hold the bracket of the clamp upwards. Now close the connector lever again so the bracket rests on the tab. Then slowly slide the Runner Handbike out of the adapter and park it safely.

NOTE!:
The coupling/uncoupling steps may differ slightly, depending on the type of stand you have chosen. This manual assumes the bicycle stand in the most basic version. See the instruction videos on the website www.roam.nl

Coupling using the Klick-Fix system
Coupling the handbike with the Klick-Fix coupling works as follows: make sure that the Klick-Fix clamps are open. Then roll your wheelchair into the Klick-Fix clamps. Apply the wheelchair brake. Close the Klick-Fix clamps one by one. Then push the stem diagonally upwards and forwards, away from you. This lifts the front wheels off the ground. Keep pressure on the stem and, with your other hand, operate the lever in the middle upwards.

Uncoupling using the Klick-Fix system
Apply the wheelchair brake. Keep pressure on the stem with one hand and then, with the other hand, operate the lever in the middle downwards. The front wheels of the wheelchair will return to the ground. Open the Klick-Fix clamps one by one. Release the wheelchair brake and roll backwards.

NOTE!:
The steps for coupling and uncoupling may differ slightly, depending on any options you have selected. This manual assumes the standard configuration in the most basic version. Please refer to the instruction videos on the website www.roam.nl.

Controls

Brakes

DANGER!
It is important that you fully understand how the braking system of your handbike works. If the braking system is not used correctly, you may lose control or it may lead to serious consequences (accident, serious injury).
The handbike is equipped with a disc brake. If, in the user's opinion, the brake does not function properly, contact the dealer immediately. The Roam dealer must resolve any defects. It is prohibited to use the handbike in any way if the braking system is not functioning.

WARNING!
1. The disc brake can become hot; contact the Roam dealer immediately if there is doubt or improper brake function.
2. Do not come into contact with the disc brake, frame parts or housing.
3. Always check your brake before riding.
4. It is prohibited to use oil or grease near the disc brake system.

Coaster brake version
Drive is provided by turning the pedals forwards. The brake is activated by turning the grips in the opposite direction: backwards!

Handbrake version
If the handbike is fitted with handbrakes on the twist grips, turning backwards as with the coaster brake does not apply; instead, simply squeeze the brake levers.

WARNING!
If the slope gradient is more than 5%, the brake can become hot and may no longer brake effectively! In this situation you must keep your speed limited (no more than 20 km per hour).

Gears

NOTE!:

The bike is supplied as standard with 7 gears or (optionally) 12 gears. The gears may only be used while riding. To shift while riding, keep the grips still.
Shift one gear at a time and reduce the force you apply to the grips while shifting.
Depending on the version chosen, the gears can be operated on the grip or on the stem.

Lighting

CAUTION!

The handbike is equipped as standard with front and rear lights. Switch on the front light by pressing the button on top of the headlamp.
The rear light must be mounted on the back of the wheelchair frame using the supplied plastic bracket with elastic mounting.
Switch on the rear light by pressing the button on the rear light unit.

Speed

NOTE ! :

The handbike is supplied as standard with 5 speed levels, which can be selected via the display. The maximum speed of the handbike is 24 km/h, unless it has been set differently at the request of your dealer or yourself.

WARNING !

- Adjust your speed to the circumstances, including the type of road surface, weather conditions and traffic situation.
- Sudden changes of direction at high speed may result in loss of control of the handbike.

Battery and charger

WARNING!

Read and follow the safety instructions before using the battery and before starting the charging process.

DANGER!

The lithium-ion battery contains chemicals that can cause dangerous reactions if the safety instructions specified in this manual are ignored.

CAUTION! Roam Special Cycles accepts no liability for damage resulting from failure to comply with these instructions.

Safety instructions and precautions for using the battery

1. Before using the Handbike for the first time, fully charge the battery.
2. The battery may only be used between -15°C and 45°C.
3. Do not expose the battery to heat/open fire. External heat may cause the battery to explode.
4. In the (unlikely) event that the battery overheats or catches fire, the battery must absolutely NOT come into contact with water or other liquids.
5. Do not open or dismantle the battery. Unauthorised opening or intentional damage to the battery poses a risk of serious injury. In addition, any warranty claim becomes invalid if the battery is opened.
6. Never connect the battery terminals in the charging socket with metal objects and ensure that the terminals never come into contact with metal objects.
7. The battery's service life depends, among other things, on where it is stored. Never leave the battery for long periods in a hot location. In particular, the boot/trunk of a car parked in the sun can be dangerous. Use the boot/trunk for transport only and not as storage space for the battery.
8. The battery must not be subjected to mechanical impact. If this happens, the battery must be inspected by the manufacturer. Contact your dealer about this. Do not continue riding with a damaged battery.
9. Always ensure that the battery is stored dry and clean.
10. When charging, place the battery on a fire- and heat-resistant surface. Do not place easily flammable objects close to the battery.
11. Do not place the charger on a wet surface or outside in the rain while charging.
12. Do not use the battery and charger if plugs or connectors are wet.
13. Improper use can lead to overheating, fire or explosion.
14. The battery must be charged to at least 20% before use.
15. During use, ensure the battery does not discharge below 20%.
16. Store the battery fully charged; if it is not used for a longer period, repeat this every 3 months.

17. The ideal storage temperature is between +10°C and +25°C and a humidity between 20% and 75%. This ensures optimal battery service life.
18. Only use original Roam chargers to charge the battery.
19. Charge the battery in a safe place under supervision; never charge at night due to the risk of fire in case of overheating.
20. Replace damaged chargers via your Roam dealer.
21. Do not store batteries in the sun or behind a window in direct sunlight.
22. When the charging cycle is complete, first unplug the charger from the wall socket before removing the battery from the charger.
23. Ensure sufficient ventilation after charging the battery.

*We recommend charging the battery in the optional case. We also recommend storing the battery in the fireproof case.

Transport
If your battery is defective, take it personally to your specialist dealer; sending lithium-ion batteries by post or other means of transport is subject to strict laws and regulations. Once again, we advise you to contact your dealer.
As transport regulations change every year, we strongly recommend consulting your travel agency, airline or shipping company before travelling, so that you know the current applicable regulations. A defective battery must not be taken on an aircraft or placed in your luggage.

REMARQUE!: Discuss transport with your dealer prior to shipment.

Using the Roam Smart charger

The MDA BC257360020 is a Li-ion battery charger with:
- Input: 100–240V AC
- Output: approx. 42V DC, 2A
- Application: charging 36V lithium-ion batteries

Use (step-by-step)
1. Check the battery
 - Only charge if it is compatible with the charger (Li-ion).
 - Do not use a damaged battery.
2. Connect the battery first
 - Connect the charger plug to the battery.
3. Then plug the charger into the mains socket
4. Indicator lights
 - Red = charging in progress
 - Green = fully charged (may vary slightly depending on the version)
5. After charging
 - First unplug from the mains socket
 - Then disconnect from the battery

Important safety instructions

- For indoor use only (no rain or moisture)
- Do not leave the charger unattended during use
- Do not open or attempt to repair yourself
- Use only with compatible batteries
- Keep out of reach of children
- Do not use in case of:
 - smoke
 - unusual smell
 - damaged cable
- Unplug from the mains socket before cleaning

Protective functions (automatic)

This charger typically has built-in protections such as:

- Overcharge protection
- Short-circuit protection
- Overheating protection
- Overvoltage protection

Tips for optimal use

- Fully charge the battery before first use
- Use the battery regularly (prevents loss of capacity)
- Do not store batteries empty
- Preferably charge at room temperature

Common problems

The charger does not work

- Check the mains socket
- Check the cables

The battery does not charge

- Possibly incorrect battery
- Defective or deeply discharged battery

The charger becomes warm

- Slightly warm = normal
- Very hot = stop immediately

Battery 36V/14.4Ah/518.4Wh

This manual describes:

- Product features and performance
- Correct use of the battery
- Safety instructions for lithium-ion battery packs

Product overview

Main parts

The battery consists of the following components:

- Lock incl. 3 keys
- LED indicator
- On/off button
- Discharge connector (to the bike)
- Locking lever (drop protection)
- Charging connector
- Warranty seal
- Product label

Technical specifications

- Battery type: Lithium-ion
- Nominal voltage: 36V
- Capacity: 14.4Ah

Use

Switching the battery on and off

- Press the on/off button to activate the battery
- The LED indicators show the charge level

Installing the battery

- Place the battery in the bike's holder
- Ensure it clicks firmly into place

Charging

First connect the charger to the battery.

1. Then plug the charger into the wall socket
2. Charging starts automatically
3. After fully charging:
 - First unplug the charger from the wall socket
 - Then disconnect the battery

LED indicator

- Multiple LEDs lit = charge level
- Flashing LEDs = charging in progress
- All LEDs lit = fully charged
 (The exact pattern may vary slightly depending on the version.)

Safety instructions

General safety

- Use only the specified charger
- Do not disassemble or modify the battery
- Protect against water and fire

- Avoid impacts, drops and damage

Safe charging
- Charge indoors only, in a dry place
- Do not cover the battery while charging
- Stop using immediately if you notice:
 - unusual smell
 - overheating
 - deformation

Safe use
- Prevent short circuits
- Do not expose to temperatures above 60°C
- Keep out of reach of children

Storage
- Store in a cool, dry place
- Ideal state of charge: 40–60%
- Recharge at least every 3 months
- Do not store empty

Maintenance
- Keep connectors clean and dry
- Do not use aggressive cleaning agents
- Check regularly for damage

Common problems
Battery does not charge
- Check the charger and connection
- Check the wall socket
- The battery may be deeply discharged

Battery does not work on the bike
- Check that it is switched on
- Check that it is installed correctly

Battery becomes too hot
- Stop using it immediately
- Allow the battery to cool down

Disposal
- Do not dispose of with household waste
- Follow local regulations for lithium batteries
- Take it to an authorised collection point

Display

CAUTION!

Do not connect or disconnect the display while it is energised. This may damage electronic components.

Due to product updates, the display may differ slightly from this manual, without affecting normal use.

Specifications:

Voltage: 12V–60V
Temperature: -20°C to +50°C
Waterproof rating: IP67
Weight: 67 g
Protocol: CAN-bus
Bluetooth: 4.0 + BLE
USB: Type-C
Certification: EC / ROHS / REACH (customisable)

Main interface

1. Headlight status (on/off)
2. Navigation (via app)
3. Maintenance reminder
4. Cruise control indicator
5. Bluetooth status
6. Battery level
7. Current speed
8. Gear position

9. Trip distance
10. Total distance
11. Average speed

Buttons
The display has 3 buttons:
- On/Off
- “+”
- “–”

General operation
On/Off
- Press and hold = on/off
- Default password: 0000

Assistance level
- “+” or “–” = adjust level (0–5)
- 0 = no assistance applied; the motor rotates freely without any resistance
- 5 = maximum assistance
- When braking or when pedalling stops, the assistance stops
- Reversing can only be done manually by operating the wheelchair push rims
- The wheelchair and handbike can be pushed at all times; the motor does not generate any resistance. The wheel rotates freely at all times.

Lighting
- Press and hold “+” = lights on/off
- The rear light must be switched on manually. It is mounted on the backrest of the wheelchair.

Walk assist (boost)
- In level 0, press “–” once more
- Then press and hold “–” to activate
- Max. speed: 6 km/h

Cruise control
- While riding, press and hold “–”
- The bike will maintain a constant speed

Display settings

Language
Choice:
- English
- German
- Spanish

Maximum speed options:
Choice of 4 firmware types depending on the selected setting and the age of the user:

- Maximum speed 10 km/h in assistance level 5
- Maximum speed 15 km/h in assistance level 5
- Maximum speed 20 km/h in assistance level 5
- Maximum speed 24 km/h in assistance level 5

Password
Enable/disable

- Press “+” and “–” for 3 seconds
- Enter password

Change password

- Enter old password
- Enter new password
- Confirm

Please note: the default is 0000

Bluetooth lock

- On/off via the menu
- Only active when enabled

Bluetooth

- Scan the QR code via the app
- After pairing, the Bluetooth icon appears

Factory reset
Restores only:

- Brightness
- Bluetooth settings
- Language
- Theme
- Password

Possible faults:

Code	Meaning	Solution
02	Brake error	Check brakes
04–05	Throttle error	Check throttle
06–07	Voltage problem	Check battery
08–10	Motor problem	Check motor
12	Current sensor error	Check controller
13–16	Battery/controller error	Stop and check
21–24	Sensor/motor errors	Check installation


29–30	Communication error	Check cables
41–48	Battery problems	Check battery
81	Bluetooth error	Check Bluetooth

Accessories

For the handbike accessories and their installation instructions, see appendix #1 et appendix #2.

Waste disposal

Below are the materials as used in the handbike, with regard to waste disposal and recycling of the handbike and the packaging.
Specific applicable legislation regarding waste disposal and recycling must be observed when the handbike has reached the end of its service life.
Aluminium: steering tube, rim
Steel: front fork, stem, motor
Stainless steel: stand, coupling frame
Plastic: grips, chain guard, (inner) tube
Packaging: cardboard, plastic

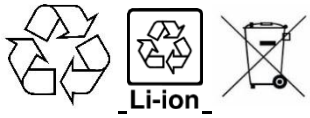Reuse

The handbike is suitable for multiple use, provided that it is cleaned and complies with the inspection procedures.
For reuse, inspect the following components for wear, damage and functionality:

- Battery: test the battery. If the measured battery capacity is above 70%, it is approved; below 70%, it is rejected. Roam can perform the battery test.
- Inspect the frame and front fork for visible cracks and paint damage.
- Inspect all brake cables and shift cables for damage and corrosion, and replace them if necessary.
- Inspect the wheel for loose spokes and play in the wheel bearings.
- Inspect the tyre for wear and replace it if necessary.
- Inspect the electrical cables for cable breaks or visible damage, and replace them if necessary.
- Inspect both brake calipers and the brake disc for proper function, including the brake pads. Replace the complete brake caliper or the brake pads in the event of poor performance or visible wear.

- Inspect the gears for proper function.
- Inspect the chain for wear using a chain wear gauge; replace it or clean it thoroughly using biodegradable degreaser.
- If present, always replace the handlebar elastic, part no. 9701052, to prevent breakage and drying out.
- Adjust the Roam adapter where necessary so that it fits the other wheelchair. NOTE: This must always be carried out by your mobility equipment dealer.
- Inspect the connector for preload and adjust it if necessary; see the mounting instruction sheet below for the connector and/or connector lever. If any parts show visible wear, always replace the complete connector lever assembly, part no. 100318, or the complete connector, part no. 9702022.
- If a Klick Fix mounting is applicable, inspect the clamping on the front frame of the wheelchair and adjust it where necessary.


Instruction sheet


1. Loosen the M6 bolts and slide them over the coupling tube; if necessary, the connector lever can be removed. Slide to the correct position and mount straight to the wheelchair.
 Tightening torque M6: 9-11Nm

2. Adjust the connector lever by adjusting the M6 nuts.

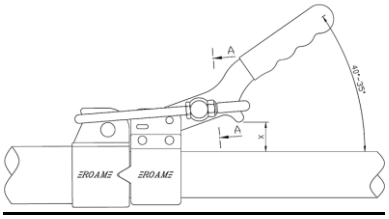3. The connector is correctly adjusted if:
 • the lever starts to clamp from 40 - 35º relative to the coupling tube and/or
 • the distance X = between 20 and 15mm

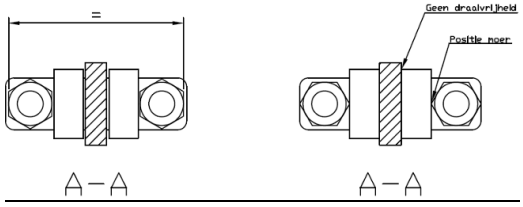4. [section A-A] Ensure that the M6 nuts are aligned with the flat side. This ensures rotational freedom at the plastic bushing and allows the lever to move smoothly.

Maintenance and cleaning

CAUTION!

Vital components such as cables, spokes and brakes are given the opportunity to settle during the initial running-in period. This period lasts approximately 1 month from the start of use. During this month, it may be necessary to make minor adjustments to the gears or brakes to ensure optimal operation. To keep the handbike in good condition, an annual service (at least once per year) is mandatory; for intensive use (more than 5,000 km per year), servicing is required twice per year. Your Roam dealer has the correct tools and the experience to adjust and lubricate concealed rotating parts such as the headset and hubs. Proper maintenance is also extremely important with regard to the warranty. Please schedule a service appointment with your Roam dealer.

Cleaning

1. Clean the handbike with a dry cloth.
2. If very dirty, clean the handbike with a soft, damp sponge and dry afterwards.
3. Do not spray the Runner Handbike with a high-pressure washer.
4. Do not use chemical cleaning agents (chlorine/bleach/soda) when cleaning the handbike.
5. Preferably store the handbike dry or at least under cover to prevent oxidation.

Tools intended for Roam dealer

The mechanic at your Roam dealer has the correct tools to maintain/repair the handbike. Original Roam parts must be used for replacement. Your Roam dealer can order these from Roam in Hedel. Ask your Roam dealer for advice.

NOTE!:

Check regularly before using the handbike:

1. Check brake function.
2. Check wear of the brake pads.
3. Check the tightening torque of all bolts and nuts.
4. Check the condition of the tyre for bumps, cracks and tread depth.
5. Check tyre pressure (recommended 1.8–2.5 bar).
6. Check the function of the front and rear lights.
7. Check the sidewall reflectors of the outer tyre.
8. Check that the grips are mounted correctly and still provide the correct grip.
9. Check the operation of the steering corrector (if fitted).
10. Check the clamping force of the connector.
11. Check the mounting of the adapter under the wheelchair frame.

WARNING!

If you have any doubts about the operation of the handbike, contact your Roam dealer immediately.

Technical specifications and features

Coupling frame:	stainless steel / aluminium
Klick-Fix coupling	stainless steel / steel
Multi stem:	stainless steel / steel
Headset:	Roam headset with lock function
Cable duct/base plate	stainless steel
Rim:	Alloy rim – 16 inch
Spokes:	Sapim
Tyre:	Schwalbe Marathon 47-305 16x1.75
Gripshifter:	Sram 3.0 comp
Rear derailleur:	Sram X4 7/9 speed
Chain:	KMC X8
Brake system:	Reverse brake Elvedes brake caliper
Front light:	LED Roam headlight (StVZO approved)
Rear light:	LED battery (StVZO approved)
Motor:	36V – 250W
Battery:	14.4Ah / 520Wh / 36V Lithium Ion
Overall width:	48 cm
Overall height:	136 cm
Overall length:	90 cm
Weight (kg):	18 kg excl. battery
Maximum load:	170 kg / 90 kg. (Klick-Fix version)
Wheel size:	16 inch
Maximum speed:	25 km/h

Tightening torques
Below is an overview of recommended tightening torques. We recommend tightening the stated parts with a torque wrench based on the specified values. Tightening parts more than specified may cause damage to your bike. For this procedure, it is best to ask your Roam dealer for help and advice.

- Expander bolt with angled wedge 12–15 Nm
- Hub axle nuts, front wheel 20–27 Nm
- Connector bolts 9–11 Nm
- Coupling frame/head tube alu clamp 38 mm 9–11 Nm
- Plastic clamp adapter (Hydac) bolts 9–11 Nm

NOTE!:
When loosening the M6 x 20 mm socket head bolts of the Roam connector on the coupling frame after making an adjustment, it is recommended to use new M6 x 20 mm bolts and refit them with assembly grease, tightening with a torque wrench to 9–11 Nm

Adjustability

The following components can be adjusted to change the user's seating position for handbike use or to change the handbike adjustment in combination with the wheelchair.

Component:	Adjustment direction:	Effect:
Coupling frame	Angular position	Coupling moment
Coupling frame on steerer tube	Height	Coupling moment
Multi stem (height)	Height	Stem higher-lower
Multi stem (length)	Longitudinal direction	Bottom bracket farther away or closer
Connector position	Longitudinal direction	Distance between handbike and wheelchair
Connector preload	Connector clamping tension	Can be adjusted tighter or looser

Possible problems

- Handbike is not aligned straight in front of the wheelchair: Loosen the 2x M6 bolts of the connector (cast part) on the coupling frame and couple the handbike to the wheelchair adapter. Align the handbike straight and then retighten these bolts alternately.
- Check whether the adapter is mounted in the centre of the wheelchair and whether the bolts are tightened to the correct torque.
- Squeaking or scraping noise at the wheel: A dragging brake caliper may be the cause; readjust the brake caliper and brake cable.
- Brake pads are worn; in that case, replace them immediately.
- Play in the connector: Readjust the connector lever. Ask your dealer to carry this out.
- Handbike does not couple properly in combination with the adapter:
- Check whether the adapter is mounted correctly under the wheelchair. It may have shifted, preventing proper coupling/uncoupling.
- The wheelchair settings may have been changed, causing the handbike to no longer couple/uncouple properly.

Warranty Terms and Conditions

Roam Special Cycles BV provides a lifetime warranty against manufacturing, material and construction defects for the frame and front fork, calculated from the date of purchase.

All parts and components manufactured and marketed by Roam Special Cycles BV are covered by a 2-year warranty against material and manufacturing defects, calculated from the date of purchase by the end user, unless stated otherwise.

The warranty excludes:
• normal wear and tear;
• tyres;
• wear parts such as brake pads, brake discs, chains, sprockets, bearings and other parts subject to normal use;
• damage resulting from improper use, insufficient maintenance, accidents, modifications to the product or use contrary to the user manual.

If a defect covered by the warranty occurs within the applicable warranty period, Roam Special Cycles BV will, at its own discretion, repair or replace the relevant part. Any shipping or transport costs for warranty parts will not be charged, unless parts unrelated to the warranty claim are shipped in the same consignment.


Warranty claims from end users must be submitted through an official Roam dealer. The relevant Roam dealer registers the product with Roam Special Cycles BV by email, stating a clear description of the complaint, the packing slip or invoice number and the relevant part. After registration, the Roam dealer sends the relevant part, together with a fully completed RMA form, to Roam Special Cycles BV.

After receipt by Roam Special Cycles BV, the warranty claim will be assessed and processed within 3 to 5 working days, unless communicated otherwise.

Please note: a copy of the original proof of purchase or packing slip and a fully completed RMA form must always be submitted. This form can be downloaded from www.roam.nl.

The warranty applies exclusively to the first owner of the product and is not transferable, unless agreed otherwise in writing.

Frame sticker

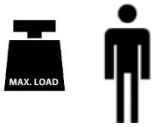TYPE:	Product name/SKU number.
	Maximum user weight.
CE	CE conformity marking.
MD	Medical device.
UDI	Unique device identification.
SN	Serial number.
	Name and address of the manufacturer.
 YYYY-MM-DD	Production date.
	Consult the instructions for use.

Manufacturer contact details

Roam Special Cycles B.V.
Buitenkamp 20
5321 KB Hedel
T: 073-5997042
E-mail: info@roam.nl
W: www.roam.nl

Closing statement

Thank you for choosing a product from Roam Special Cycles. This product has been developed in accordance with high quality and safety requirements. With correct use and maintenance, it will contribute to reliable and comfortable mobility.
For additional information or support, please first contact your supplier/dealer or Roam Special Cycles.

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