Handleiding
Directions for use
Notice d'utilisation
Gebrauchsanweisung

E-bike 20''

CE

MD

ROAM Special Cycles B.V.
Buitenkamp 20
5321KB Hedel
The Netherlands
T: 0031 (0)73-5997042
E-mail: info@roam.nl
W: www.roam.nl

Inleiding/Introduction/Introduction/Einleitung

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Table of Contents

Introduction	2
Table of Contents	34
Definitions	35
Delivery inspection	35-36
General	36
Applications	37
Regulations	38
Use	38-42
Coupling and uncoupling the E-bike 20''	42-43
Controls	43-44
Battery and charger	45-50
Display	50-53
Accessories	53
Waste disposal	53
Reuse	54-55
Maintenance and cleaning	56
Technical specifications and features	57-58
Possible problems	59
Warranty terms and conditions	59-60
Frame sticker	60
Manufacturer contact details	61
Closing remarks	61
Appendices	119-131

Definitions

Always observe the definitions below!

Word	Definition
 **DANGER!**	The user is alerted to the risk of serious injury or death if the advice is not followed.
 **WARNING!**	The user is alerted to the risk of injury if the advice is not followed.
 **CAUTION!**	The user is alerted to the risk of damage to the equipment if the advice is not followed.
NOTE:	General advice or best practice
	Reference to additional documentation

Delivery check

Check whether the delivery matches the delivery note. If transport damage is found upon delivery, report this (within 24 hours) to your Roam dealer, supported with photos. Roam trusts that the E-bike 20'' meets your expectations and wishes you enjoyable ride.

Scope of Delivery

- E-bike 20''
- Roam Adapter (if applicable)
- Charger
- Battery (check the battery type on the packing slip and on the label of the battery housing)
- Separate rear light
- Manual
- Other components depend on the options selected.

NOTE!:
Before delivery of your Roam, all checkpoints have been inspected using a checklist. Certain parts must be readjusted after some time. We recommend having this carried out by your Roam dealer.

NOTE!:
Write down the details of the dealer who supplied your product here. In the event of malfunctions, contact your dealer.

Dealer name:

Dealer phone number:

Dealer email address:

Date:

General

NOTE!:
Failure to follow these instructions may result in physical injury, damage to the product, or damage to the surroundings! A notice to the user and/or patient: If this product has been involved in a serious incident, it must be reported to Roam Special Cycles B.V. and the competent authority of the Member State in which the user and/or patient is established.
The expected service life of the product is 7 years.
It may be possible to combine this medical device with one or more other medical devices or another product. For more information, please contact **Roam Special Cycles B.V.**

No modifications to the E-bike 20''

WARNING!
Roam cannot be held liable for incorrect use or incorrect information regarding the use of the E-bike 20'', and any resulting damage, unsafe situations/actions, accidents or incidents.
No modifications may be made to the E-bike 20''.
Roam is not responsible for damage to the E-bike 20'' caused by modifications.
The speed of our products must not be increased or altered in any way. If the product is modified, the producer/manufacturer is in no way responsible or liable.

Applications

Indications
The E-bike 20'' provides relief and compensation in cases of physical immobility, disability, or another limitation. The E-bike 20'' is a Class I medical device, in accordance with Annex VIII of the MDR. Suitable for both adults and children with a body height from 1.10 metres to 2.10 metres and an age of 5 years and older.

Contra-indications
The E-bike 20'' must not be used without supervision in cases of:

- Perceptual disorders
- Loss of both arms
- Inability to adequately control braking and steering
- Recent surgery to the arms, shoulders and hands
- Severe fatigue or limited capacity

Additional
NOTE!:
The user must be able to operate the E-bike 20'' properly.
The user of the E-bike 20'' must be in good mental health and must not be physically restricted in a seated position.
The rider's hearing and eyesight must be adequate for using the E-bike 20''.
The user must read and understand the manual thoroughly.
If the user is so limited that they have vision/hearing problems and/or cannot read and understand the manual, parents/carers/supervisors are responsible for applying the manual and for the end user's use of the E-bike 20''.
The E-bike 20'' is for personal use. Only the wheelchair user is permitted to use the E-bike 20''. The user is responsible and liable for correct use as stated in this manual.

Additional information for the dealer
NOTE!:
It is prohibited to use the E-bike 20'' under the influence of alcohol/drugs or mind-altering medication.
Read the instructions for use of the E-bike 20''. Instruct the user on the safe use of the E-bike 20'' and any additionally fitted components. Inform users about specific warnings that must be read, understood and observed. When delivery/purchase of a E-bike 20'' is being considered, take into account the following: body size, weight, physical and mental condition, age, living conditions and environment.

Important: Do not use your E-bike 20'' until you have read and understood these instructions for use.

Regulations

NOTE!:

The user must at all times comply with the laws and regulations applicable in the country of use regarding road safety and the use of this product. The manufacturer accepts no liability for damage resulting from failure to comply with these regulations.

NOTE!:

For users under 16 years of age, it is recommended that the E-bike 20'' is used under the supervision of a parent or carer.
The user may use: the footpath, the (moped) cycle path and the roadway.
It is prohibited to use motorways or expressways.
Comply with the traffic rules applicable to the road you are using. On the pavement, the maximum speed is 6 km/h.
Electronic devices such as mobile phones, navigation systems, tablets or music players must not be held while riding; Roam strongly discourages using them while riding. The same rules apply as for texting, calling or listening to music while cycling.

Use

WARNING!

1. Adjust your speed to the conditions (e.g. heavy traffic, poor road surface, obstacles and roads alongside water). Take speed bumps and other uneven surfaces at an appropriate speed.
2. Keep both hands on the handlebar/hand grips while riding and do not allow yourself to be distracted. Never put your hands into moving parts.
3. Be extra alert near obstacles, raised surfaces, railway crossings, tram rails, pedestrian crossings and when riding along busy roads or water.
4. Use lights at dusk, in the dark or in poor weather conditions. Roam also recommends using lights during the day.
5. Roam recommends that wheelchair users wear a helmet when using the E-bike 20''.
6. Make sure you are well prepared while on the road: keep the contact details of your Roam dealer at hand (e.g. in your mobile phone).
7. Protect yourself during longer rides against cold and precipitation in winter and against excessive heat in summer.
8. Before every ride, check that the adapter and/or clamping systems under the wheelchair and/or on the E-bike 20'' are correctly secured. For removable adapters, check the star knobs/quick-release fasteners.
9. Have maintenance, adjustments, replacement of parts and repairs carried out by your Roam dealer.
10. Take right-angle and sharp turns at an appropriate speed to prevent the risk of tipping.

Slopes and brakes

WARNING!


1. In some cases, if this manual is not followed correctly, a dangerous situation may arise for the E-bike 20'' user when riding on slopes.
2. We recommend not using the E-bike 20'' on slopes of 5% or more. On very long or steep slopes, the brakes may overheat and the motor may be overloaded, which can cause defects.
3. The point above regarding slopes also applies to descents.
4. Counterweights (optional) can be fitted, generating more wheel load and therefore more grip on the front wheel. This advantage of increased wheel load is beneficial for uphill slopes, not for descents (speed may increase).
5. Be aware that speed increases when descending a slope.
6. Please note: On roads with oil, mud, sand or leaves on the surface, the braking distance may be longer. This also applies to uneven road surfaces (potholes/holes) and winter conditions.

Front wheel traction loss

DANGER!

When using the E-bike 20'', the front wheel may slip under certain conditions, particularly when riding on slopes or on slippery surfaces. This phenomenon can affect riding performance and steering control. Although the risk of slipping can be reduced through correct use and adjustment, it cannot be completely eliminated.

- Ensure that the tyre is not overinflated (2 to 2.5 bar)
- Ride in the highest assistance mode (mode 5)
- Maintain momentum. Build up speed in advance, not only once you are already on the slope
- Where possible, shift your body slightly forward
- Install counterweights; these reduce slipping
- If necessary, fit an additional set of counterweight plates to the counterweights
- If necessary, install a cargo carrier so that carried items are mounted on the bike and not on the back of the wheelchair
- Ensure that the lift of the wheelchair's front wheels is no more than 5 cm
- Use a more active E-bike 20'' adjustment, with the unit positioned further forward
- Install a second rear axle on the wheelchair (not possible on every wheelchair; consult your dealer/wheelchair manufacturer)
- Keep your speed up until you reach the top, and only then slow down.
- Take extra care on wet, loose or slippery surfaces

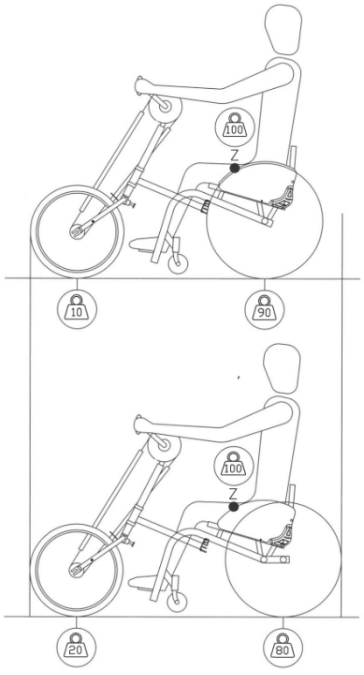*The image shows an example of a second rear axle; the pressure on the front wheel increases.

Corners

**DANGER!**

Right-angle and sharp turns must be taken at an appropriate speed, depending on the conditions (road surface and weather) and physical capabilities, to prevent the risk of tipping. Make sure that your upper legs or upper body are not struck by the handles, as this may restrict the steering angle.

Maximum user weight

**DANGER!**

Never exceed the maximum user weight of 170 kg. This is separate from the maximum user weight of the wheelchair; please consult the relevant wheelchair supplier/manufacturer. Exceeding the maximum load may cause damage to the E-bike 20'', or you may fall over or tip, lose control, and cause serious injury to the user and other people.

Choking hazard

**WARNING!**

This mobility aid contains small parts which, under certain circumstances, may pose a choking hazard for small children. The product shown and described in this manual may not be exactly the same in every detail as your model. However, all instructions are important regardless of differences in detail. The manufacturer reserves the right to make changes without prior notice to information regarding weight, dimensions or other technical data as stated in this manual.

Wheelchair

**WARNING!**

If the wheelchair balance point is adjusted to be sharp or tip-sensitive, the E-bike 20'' will have reduced wheel load and the braking distance will be much longer. Reduced wheel load can cause a risk of slipping or wheel spin, including when riding uphill.

NOTE:

The Roam adapter increases the weight of the wheelchair by 1 to 2 kg. This depends on the type of wheelchair and the associated fastening components of the Roam adapter. This has little to no effect on the driving characteristics of the wheelchair itself when the E-bike 20'' is not in use.

If the user wishes to remove the adapter from the wheelchair, this can be done by loosening the M6 star knobs (counterclockwise) and removing the adapter from the wheelchair with a forward movement. This applies to a removable adapter.
In the Klick-Fix (front frame mounting) version, nothing additional is attached to the wheelchair for coupling.

The wheelchair can be fully restored to its original condition after removal of the Roam E-bike 20'' and the associated adapter construction. In doing so, all fastening materials and mounting components must be completely removed from the wheelchair frame.
In a version with the Klick-Fix system, no permanent components or fasteners are mounted on the original wheelchair frame. As a result, the original configuration of the wheelchair remains unchanged.

Coupling and uncoupling the E-bike 20''

Coupling using the adapter system
Coupling the E-bike 20'' to your wheelchair is very easy. First apply the wheelchair brake; then close the connector lever on the E-bike 20'' (downwards) and slide the E-bike 20'' with the coupling frame as far as possible into the adapter on your wheelchair.
Now open the connector lever on the E-bike 20'' (upwards) so that the bracket falls over the tab of the adapter; then close the lever again. The E-bike 20'' is now coupled. Release the wheelchair brake and you are ready to ride.

NOTE!:
The coupling/uncoupling steps may differ slightly, depending on the type of stand you have chosen. This manual assumes the bicycle stand in the most basic version. See the instruction videos on the website www.roam.nl

Uncoupling using the adapter system
Uncoupling the E-bike 20'' is also very easy. First apply the wheelchair brake. Then hold the E-bike 20'' (preferably by holding the handlebars); with your other hand open the connector and hold the bracket of the clamp upwards. Now close the connector lever again so the bracket rests on the tab. Then slowly slide the E-bike 20'' out of the adapter and park it safely.

NOTE!:
The coupling/uncoupling steps may differ slightly, depending on the type of stand you have chosen. This manual assumes the bicycle stand in the most basic version. See the instruction videos on the website www.roam.nl

Coupling using the Klick-Fix system
Coupling the E-bike 20'' with the Klick-Fix coupling works as follows: make sure that the Klick-Fix clamps are open. Then roll your wheelchair into the Klick-Fix clamps. Apply the wheelchair brake. Close the Klick-Fix clamps one by one. Then push the stem diagonally upwards and forwards, away from you. This lifts the front wheels off the ground. Keep pressure on the stem and, with your other hand, operate the lever in the middle upwards.

Uncoupling using the Klick-Fix system
Apply the wheelchair brake. Keep pressure on the stem with one hand and then, with the other hand, operate the lever in the middle downwards. The front wheels of the wheelchair will return to the ground. Open the Klick-Fix clamps one by one. Release the wheelchair brake and roll backwards.

NOTE!:
The steps for coupling and uncoupling may differ slightly, depending on any options you have selected. This manual assumes the standard configuration in the most basic version. Please refer to the instruction videos on the website www.roam.nl.

Controls

Brakes

**DANGER!**
It is important that you fully understand the operation of the braking system of your E-bike 20''. If the braking system is not used correctly, this may cause you to lose control of the handlebars or may lead to serious consequences (accident, serious injuries).
The E-bike 20'' is equipped with a disc brake. If, in the user's opinion, the brakes do not function properly, the dealer must be contacted immediately. The Roam dealer must resolve any defects. It is prohibited to use the E-bike 20'' if the braking system is not functioning in any form.

**WARNING!**

1. The disc brake may become hot; contact the Roam dealer immediately if there is any doubt or improper operation when the brake is functioning.
2. Do not come into contact with the disc brake, structure or housing.
3. Always check your brake before riding.
4. It is prohibited to use oil or grease near the disc brake system.

Braking system

The E-bike 20'' is equipped with a hand brake on both sides. Each brake lever has its own brake caliper operating on the same brake disc.

Easy E-bike tilt handlebar (optional)

If the E-bike 20'' is equipped with the Easy E-bike tilt handlebar, the brake is activated when you push the handlebar away from you. Throttle control is activated by pulling the handlebar towards you. The extent to which braking and throttle control are activated (hard braking/high throttle) depends on how far the handlebar is moved forwards or backwards. This handlebar is optional on the E-bike 20''.

**WARNING!**
If the slope gradient is more than 5%, the brake can become hot and may no longer brake effectively! In this situation you must keep your speed limited (no more than 20 km per hour).

Lighting

CAUTION!

The E-bike 20'' is equipped as standard with front and rear lighting. Switch on the front light by pressing the button on top of the headlight.
The rear light must be mounted on the back of the relevant wheelchair frame using the supplied plastic bracket with elastic mounting.
Switch on the rear light by pressing the button on the rear light unit.

Speed

NOTE!:

The E-bike 20'' is supplied as standard with 5 speed levels, which can be selected via the display. The maximum speed of the E-bike 20'' is 20 km/h, unless this has been set differently at the request of your dealer or yourself.

CAUTION!

- Adjust your speed to the circumstances, including the type of road surface, weather conditions and traffic situation.
- Sudden changes of direction at higher speed may result in loss of control of the E-bike 20''.

Transport

When transporting the E-bike 20'', it is recommended to use a suitable box with sufficient protective material to prevent damage.

The E-bike 20'' may be taken on an aircraft; also use a suitable box and protective material for this. Always ask the airline for advice and information.

NOTE!:

Discuss transport to your dealer with your dealer before shipping.

Battery and charger

WARNING!

Read and follow the safety instructions before using the battery and before starting the charging process.

DANGER!

The lithium-ion battery contains chemicals that can cause dangerous reactions if the safety instructions specified in this manual are ignored.

CAUTION! Roam Special Cycles accepts no liability for damage resulting from failure to comply with these instructions.

Safety instructions and precautions for using the battery

1. Before using the E-bike 20'' for the first time, the battery must be fully charged.
2. The battery may only be used between -15°C and 45°C.
3. Do not expose the battery to heat/open fire. External heat may cause the battery to explode.
4. In the (unlikely) event that the battery overheats or catches fire, the battery must absolutely NOT come into contact with water or other liquids.
5. Do not open or dismantle the battery. Unauthorised opening or intentional damage to the battery poses a risk of serious injury. In addition, any warranty claim becomes invalid if the battery is opened.
6. Never connect the battery terminals in the charging socket with metal objects and ensure that the terminals never come into contact with metal objects.
7. The battery's service life depends, among other things, on where it is stored. Never leave the battery for long periods in a hot location. In particular, the boot/trunk of a car parked in the sun can be dangerous. Use the boot/trunk for transport only and not as storage space for the battery.
8. The battery must not be subjected to mechanical impact. If this happens, the battery must be inspected by the manufacturer. Contact your dealer about this. Do not continue riding with a damaged battery.
9. Always ensure that the battery is stored dry and clean.
10. When charging, place the battery on a fire- and heat-resistant surface. Do not place easily flammable objects close to the battery.
11. Do not place the charger on a wet surface or outside in the rain while charging.
12. Do not use the battery and charger if plugs or connectors are wet.
13. Improper use can lead to overheating, fire or explosion.
14. The battery must be charged to at least 20% before use.
15. During use, ensure the battery does not discharge below 20%.
16. Store the battery fully charged; if it is not used for a longer period, repeat this every 3 months.

17. The ideal storage temperature is between +10°C and +25°C and a humidity between 20% and 75%. This ensures optimal battery service life.
18. Only use original Roam chargers to charge the battery.
19. Charge the battery in a safe place under supervision; never charge at night due to the risk of fire in case of overheating.
20. Replace damaged chargers via your Roam dealer.
21. Do not store batteries in the sun or behind a window in direct sunlight.
22. When the charging cycle is complete, first unplug the charger from the wall socket before removing the battery from the charger.
23. Ensure sufficient ventilation after charging the battery.

*We recommend charging the battery in the optional case. We also recommend storing the battery in the fireproof case.

Transport
If your battery is defective, take it personally to your specialist dealer; sending lithium-ion batteries by post or other means of transport is subject to strict laws and regulations. Once again, we advise you to contact your dealer.
As transport regulations change every year, we strongly recommend consulting your travel agency, airline or shipping company before travelling, so that you know the current applicable regulations. A defective battery must not be taken on an aircraft or placed in your luggage.

REMARQUE!: Discuss transport with your dealer prior to shipment.

Using the Roam Smart charger

The MDA BC257360020 is a Li-ion battery charger with:
- Input: 100–240V AC
- Output: approx. 42V DC, 2A
- Application: charging 36V lithium-ion batteries

Use (step-by-step)
1. Check the battery
 - Only charge if it is compatible with the charger (Li-ion).
 - Do not use a damaged battery.
2. Connect the battery first
 - Connect the charger plug to the battery.
3. Then plug the charger into the mains socket
4. Indicator lights
 - Red = charging in progress
 - Green = fully charged (may vary slightly depending on the version)
5. After charging
 - First unplug from the mains socket
 - Then disconnect from the battery

Important safety instructions

- For indoor use only (no rain or moisture)
- Do not leave the charger unattended during use
- Do not open or attempt to repair yourself
- Use only with compatible batteries
- Keep out of reach of children
- Do not use in case of:
 - smoke
 - unusual smell
 - damaged cable
- Unplug from the mains socket before cleaning

Protective functions (automatic)

This charger typically has built-in protections such as:

- Overcharge protection
- Short-circuit protection
- Overheating protection
- Overvoltage protection

Tips for optimal use

- Fully charge the battery before first use
- Use the battery regularly (prevents loss of capacity)
- Do not store batteries empty
- Preferably charge at room temperature

Common problems

The charger does not work

- Check the mains socket
- Check the cables

The battery does not charge

- Possibly incorrect battery
- Defective or deeply discharged battery

The charger becomes warm

- Slightly warm = normal
- Very hot = stop immediately

Battery 36V/14.4Ah/518.4Wh

This manual describes:

- Product features and performance
- Correct use of the battery
- Safety instructions for lithium-ion battery packs

Product overview

Main parts

The battery consists of the following components:

- Lock incl. 3 keys
- LED indicator
- On/off button
- Discharge connector (to the bike)
- Locking lever (drop protection)
- Charging connector
- Warranty seal
- Product label

Technical specifications

- Battery type: Lithium-ion
- Nominal voltage: 36V
- Capacity: 14.4Ah

Use

Switching the battery on and off

- Press the on/off button to activate the battery
- The LED indicators show the charge level

Installing the battery

- Place the battery in the bike's holder
- Ensure it clicks firmly into place

Charging

First connect the charger to the battery.

1. Then plug the charger into the wall socket
2. Charging starts automatically
3. After fully charging:
 - First unplug the charger from the wall socket
 - Then disconnect the battery

LED indicator

- Multiple LEDs lit = charge level
- Flashing LEDs = charging in progress
- All LEDs lit = fully charged
 (The exact pattern may vary slightly depending on the version.)

Safety instructions

General safety

- Use only the specified charger
- Do not disassemble or modify the battery
- Protect against water and fire

- Avoid impacts, drops and damage

Safe charging
- Charge indoors only, in a dry place
- Do not cover the battery while charging
- Stop using immediately if you notice:
 - unusual smell
 - overheating
 - deformation

Safe use
- Prevent short circuits
- Do not expose to temperatures above 60°C
- Keep out of reach of children

Storage
- Store in a cool, dry place
- Ideal state of charge: 40–60%
- Recharge at least every 3 months
- Do not store empty

Maintenance
- Keep connectors clean and dry
- Do not use aggressive cleaning agents
- Check regularly for damage

Common problems
Battery does not charge
- Check the charger and connection
- Check the power socket
- The battery may be deeply discharged

Battery does not work on the bike
- Check that it is switched on
- Check that it is installed correctly

Battery becomes too hot
- Stop using it immediately
- Allow the battery to cool down

Disposal
- Do not dispose of with household waste
- Follow local regulations for lithium batteries
- Take it to an authorised collection point

Display

CAUTION!

Do not connect or disconnect the display while it is energised. This may damage electronic components.

Due to product updates, the display may differ slightly from this manual, without affecting normal use.

Specifications:

Voltage: 12V–60V
Temperature: -20°C to +50°C
Waterproof rating: IP67
Weight: 67 g
Protocol: CAN
Bluetooth: 4.0 + BLE
USB: Type-C
Certification: EC / ROHS / REACH (customisable)

Main interface

1. Headlight status (on/off)
2. Navigation (via app)
3. Maintenance reminder
4. Cruise control indicator
5. Bluetooth status
6. Battery level
7. Current speed
8. Gear position

9. Trip distance
10. Total distance
11. Average speed

Buttons
The display has 3 buttons:
- On/Off
- “+”
- “–”

General operation
On/Off
- Press and hold = on/off
- Default password: 0000

Speed level
- “+” or “–” = adjust level (0–5)
- 0 = no drive assistance; the motor rotates freely without any resistance
- 5 = maximum assistance
- When braking or when you stop turning the cranks, the assistance stops
- Reverse can be selected by pressing the “–” button until “R” appears on the display. Use the “+” button to select the desired speed level again for forward travel.
- The wheelchair and E-bike 20'' can be pushed at all times; the motor does not create resistance. The wheel rotates freely at all times.

Lighting
- Press and hold “+” = lights on/off
- The rear light must be switched on manually. It is mounted on the backrest of the wheelchair.

Cruise control
- During riding, press and hold “–”.
- Gradually return the throttle to the zero position.
- The bike maintains a constant speed.
- Applying the brake or throttle disengages the cruise control.

Display settings

Language
Choice:
- English
- German
- Spanish

Maximum speed options:
Choice of 3 firmware types depending on the selected setting and the age of the user:

- Maximum speed 10 km/h in assistance level 5
- Maximum speed 15 km/h in assistance level 5
- Maximum speed 20 km/h in assistance level 5

Password
Enable/disable

- Press “+” and “–” for 3 seconds
- Enter password

Change password

- Enter old password
- Enter new password
- Confirm

Please note: the default is 0000

Bluetooth lock

- On/off via the menu
- Only active when enabled

Bluetooth

- Scan the QR code via the app
- After pairing, the Bluetooth icon appears

Factory reset
Restores only:

- Brightness
- Bluetooth settings
- Language
- Theme
- Password

Possible faults:

Code	Meaning	Solution
02	Brake error	Check brakes
04–05	Throttle error	Check throttle
06–07	Voltage problem	Check battery
08–10	Motor problem	Check motor
12	Current sensor error	Check controller
13–16	Battery/controller error	Stop and check
21–24	Sensor/motor errors	Check installation


29–30	Communication error	Check cables
41–48	Battery problems	Check battery
81	Bluetooth error	Check Bluetooth

Accessories

For the E-bike 20'' accessories and their installation instructions, see appendix #1 , appendix #2 and appendix #3.

Waste disposal

Below are the materials as used in the E-bike 20'', with regard to waste disposal and recycling of the E-bike 20'' and the packaging.
Specific applicable legislation regarding waste disposal and recycling must be observed when the E-bike 20'' has reached the end of its service life.
Aluminium: steering tube, rim
Steel: front fork, stem, motor
Stainless steel: stand, coupling frame
Plastic: grips, chain guard, (inner) tube
Packaging: cardboard, plastic

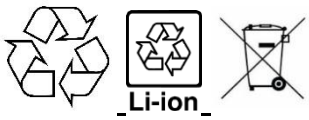Reuse

The E-bike 20'' is suitable for multiple use, provided that it is cleaned and complies with the inspection procedures.
Check the following products for reuse for wear, damage and functionality before reuse:

- Battery: test it. If the measured battery capacity is above 70%, it is approved; below 70%, it is rejected. Battery testing can be carried out by Roam.
- Check the frame and front fork for visible cracks and paint damage.
- Check the condition of all brake cables and shift cables for damage and corrosion and replace them if necessary.
- Check the wheel for loose spokes and play in the wheel bearings.
- Check the tyre for wear and replace it if necessary.
- Check the electrical cables for cable breakage or visible damage and replace them if necessary.
- Check both brake calipers and the brake disc for functionality, including the brake pads. Replace the entire brake caliper or the brake pads in the event of poor functioning or visible wear.
- Check the chain for wear with a chain wear gauge; replace it or clean it thoroughly with a biodegradable degreaser.
- If present, always replace the steering stabiliser elastic, art. 9701052, to prevent breakage and drying out.
- Adjust the Roam adapter where necessary so that it fits the other wheelchair. NOTE!: This must always be carried out by your assistive devices dealer.
- Check the connector for preload and adjust it if necessary; for this, see the instruction sheet below: connector assembly instructions and/or connector lever. If parts show visible wear, always replace the complete connector lever assembly, art. 100318, or the complete connector, art. no. 9702022.
- If a Klick Fix mounting is applicable, check the clamping on the wheelchair front frame and adjust it where necessary.


Instruction sheet


1. Loosen the M6 bolts and slide them over the coupling tube; if necessary, the connector lever can be removed. Slide to the correct position and mount straight to the wheelchair.
 Tightening torque M6: 9-11Nm

2. Adjust the connector lever by adjusting the M6 nuts.

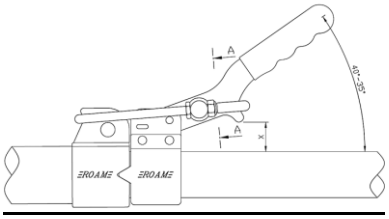3. The connector is correctly adjusted if:
 • the lever starts to clamp from 40 - 35º relative to the coupling tube and/or
 • the distance X = between 20 and 15mm

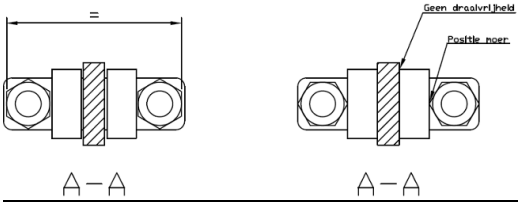4. [section A-A] Ensure that the M6 nuts are aligned with the flat side. This ensures rotational freedom at the plastic bushing and allows the lever to move smoothly.

Maintenance and cleaning

CAUTION!

Vital components such as cables, spokes and brakes are given the opportunity to settle during the initial running-in period. This period lasts approximately 1 month from the start of use. During this month, it may be necessary to make minor adjustments to the gears or brakes to ensure optimal operation. To keep the E-bike 20'' in good condition, an annual service (at least once per year) is mandatory; for intensive use (more than 5,000 km per year), servicing is required twice per year. Your Roam dealer has the correct tools and the experience to adjust and lubricate concealed rotating parts such as the headset and hubs. Proper maintenance is also extremely important with regard to the warranty. Please schedule a service appointment with your Roam dealer.

Cleaning

1. Clean the E-bike 20'' with a dry cloth.
2. If very dirty, clean the E-bike 20'' with a soft, damp sponge and dry afterwards.
3. Do not spray the E-bike 20'' with a high-pressure washer.
4. Do not use chemical cleaning agents (chlorine/bleach/soda) when cleaning the E-bike 20''.
5. Preferably store the E-bike 20'' dry or at least under cover to prevent oxidation.

Tools intended for the Roam dealer

The technician at your Roam dealer has the correct tools to maintain/repair the E-bike 20''. Original Roam parts must be used for replacement. Your Roam dealer can order these from Roam Special Cycles BV. Ask your Roam dealer for advice on this.

NOTE!:

Check regularly before using the E-bike 20'':

1. Check brake function.
2. Check wear of the brake pads.
3. Check the tightening torque of all bolts and nuts.
4. Check the condition of the tyre for bumps, cracks and tread depth.
5. Check tyre pressure (recommended 2.5–3.5 bar).
6. Check the function of the front and rear lights.
7. Check the sidewall reflectors of the outer tyre.
8. Check that the grips are mounted correctly and still provide the correct grip.
9. Check the operation of the steering corrector (if fitted).
10. Check the clamping force of the connector.
11. Check the mounting of the adapter under the wheelchair frame.

WARNING!

If you have any doubts about the operation of the E-bike 20'', contact your Roam dealer immediately.

Technical specifications and features

Adjustable coupling frame:	stainless steel / aluminium
Klick-Fix coupling:	stainless steel / steel
Adjustable stem:	Roam aluminium
Handlebar:	Roam aluminium 56 cm wide
Headset:	Roam headset with locking function
Cable duct/base plate:	stainless steel
Rim:	Ryde Andra 29 - 20 inch
Spokes:	Sapim
Tyre:	Schwalbe Marathon 47-406 20x1.75
Braking system:	Hand brakes L/R Elvedes brake calipers
Front lighting:	LED Roam headlamp (StVZO approved)
Rear lighting:	LED battery light (StVZO approved)
Motor:	36V – 250W
Battery:	14.4Ah / 520Wh / 36V Lithium-ion
Overall width:	58 cm
Overall height:	140 cm
Overall length:	90 cm
Weight in kg:	14 kg excl. battery
Maximum load:	170 kg / 90 kg (Klick-Fix version)
Wheel size:	20 inch
Maximum speed:	20 km/h

Tightening torques
Below is an overview of recommended tightening torques. We recommend tightening the stated parts with a torque wrench based on the specified values. Tightening parts more than specified may cause damage to your bike. For this procedure, it is best to ask your Roam dealer for help and advice.

- Expander bolt with angled wedge 12–15 Nm
- Hub axle nuts, front wheel 20–27 Nm
- Connector bolts 9–11 Nm
- Coupling frame/head tube alu clamp 38 mm 9–11 Nm
- Plastic clamp adapter (Hydac) bolts 9–11 Nm

NOTE!:
When loosening the M6 x 20 mm socket head bolts of the Roam connector on the coupling frame after making an adjustment, it is recommended to use new M6 x 20 mm bolts and refit them with assembly grease, tightening with a torque wrench to 9–11 Nm.

Adjustability

The following components can be adjusted to change the user's seating position for E-bike 20'' use or to change the E-bike 20'' adjustment in combination with the wheelchair.

Component:	Adjustment direction:	Effect:
Coupling frame	Angle adjustment	Coupling moment
Coupling frame on steering tube	Height adjustment	Coupling moment
Stem (height)	Height adjustment	Stem up/down
Stem (angle adjustment)	Stem steering angle	Handlebar closer/further away
Connector position	Longitudinal adjustment	Distance between the E-bike 20'' and the wheelchair
Connector preload	Connector clamping tension	Adjustable tighter or looser

Possible problems

- E-bike 20'' is not aligned straight in front of the wheelchair: Loosen the 2x M6 bolts of the connector (cast part) on the coupling frame and couple the E-bike 20'' to the wheelchair adapter. Align the E-bike 20'' straight and then retighten these bolts alternately.
- Check whether the adapter is mounted in the centre of the wheelchair and whether the bolts are tightened to the correct torque.
- Squeaking or scraping noise at the wheel: A dragging brake caliper may be the cause; readjust the brake caliper and brake cable.
- Brake pads are worn; in that case, replace them immediately.
- Play in the connector: Readjust the connector lever. Ask your dealer to carry this out.
- E-bike 20'' does not couple properly in combination with the adapter:
- Check whether the adapter is mounted correctly under the wheelchair. It may have shifted, preventing proper coupling/uncoupling.
- The wheelchair settings may have been changed, causing the E-bike 20'' to no longer couple/uncouple properly.

Warranty Terms and Conditions

Roam Special Cycles BV provides a lifetime warranty against manufacturing, material and construction defects for the frame and front fork, calculated from the date of purchase.

All parts and components manufactured and marketed by Roam Special Cycles BV are covered by a 2-year warranty against material and manufacturing defects, calculated from the date of purchase by the end user, unless stated otherwise.

The warranty excludes:
• normal wear and tear;
• tyres;
• wear parts such as brake pads, brake discs, chains, sprockets, bearings and other parts subject to normal use;
• damage resulting from improper use, insufficient maintenance, accidents, modifications to the product or use contrary to the user manual.

If a defect covered by the warranty occurs within the applicable warranty period, Roam Special Cycles BV will, at its own discretion, repair or replace the relevant part. Any shipping or transport costs for warranty parts will not be charged, unless parts unrelated to the warranty claim are shipped in the same consignment.


Warranty claims from end users must be submitted through an official Roam dealer. The relevant Roam dealer registers the product with Roam Special Cycles BV by email, stating a clear description of the complaint, the packing slip or invoice number and the relevant part. After registration, the Roam dealer sends the relevant part, together with a fully completed RMA form, to Roam Special Cycles BV.

After receipt by Roam Special Cycles BV, the warranty claim will be assessed and processed within 3 to 5 working days, unless communicated otherwise.

Please note: a copy of the original proof of purchase or packing slip and a fully completed RMA form must always be submitted. This form can be downloaded from www.roam.nl.

The warranty applies exclusively to the first owner of the product and is not transferable, unless agreed otherwise in writing.

Frame sticker

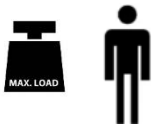TYPE:	Product name/SKU number.
	Maximum user weight.
	CE conformity marking.
	Medical device.
	Unique device identification.
	Serial number.
	Name and address of the manufacturer.
 YYYY-MM-DD	Production date.
	Consult the instructions for use.

Manufacturer contact details

Roam Special Cycles B.V.
Buitenkamp 20
5321 KB Hedel
T: 073-5997042
E-mail: info@roam.nl
W: www.roam.nl

Closing statement

Thank you for choosing a product from Roam Special Cycles. This product has been developed in accordance with high quality and safety requirements. With correct use and maintenance, it will contribute to reliable and comfortable mobility.
For additional information or support, please first contact your supplier/dealer or Roam Special Cycles.

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